


NIDEK

OPHTHALMIC YAG LASER SYSTEM YC-200

OPERATOR'S MANUAL

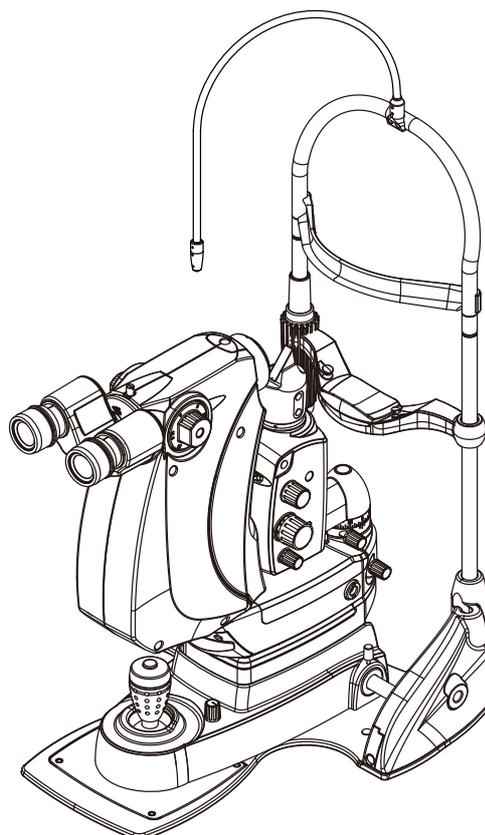

Be sure to read the SOFTWARE LICENSE AGREEMENT (page 2) before using this product.

Original instructions

NIDEK CO., LTD.

NIDEK CO., LTD. (Manufacturer)	: 34-14 Maehama, Hiroishi-cho, Gamagori, Aichi 443-0038, JAPAN Telephone: +81-533-67-6611 URL: https://www.nidek.com/
NIDEK INC. (United States Agent)	: 2040 Corporate Court, San Jose, CA 95131, U.S.A. Telephone: +1-800-223-9044 (USA Only) URL: https://usa.nidek.com/
NIDEK S.A. (EU Authorized Representative)	: Ecoparc, rue Benjamin Franklin, 94370 Sucy En Brie, FRANCE



2024-04-12
YC020-P902-B2
Printed in Japan
© 2018 NIDEK CO., LTD.

Before Use

This operator's manual describes the operating procedure, safety precautions, maintenance, and specifications of the NIDEK Ophthalmic YAG Laser System YC-200. Be sure to read the operator's manual prior to use of the device to understand the safety precautions and operating procedures thoroughly.

Keep this manual handy for reference.

If you encounter any problems or have questions about the device, please contact NIDEK or your authorized distributor.

IMPORTANT - READ CAREFULLY

THIS AGREEMENT APPLIES TO THE NIDEK SOFTWARE AND ACCOMPANYING DOCUMENTS. PLEASE READ THIS AGREEMENT CAREFULLY AND THOROUGHLY BEFORE USING SOFTWARE.

SOFTWARE LICENSE AGREEMENT

This SOFTWARE LICENSE AGREEMENT (this "Agreement") is an agreement between you, whether person or legal entity, and NIDEK CO., LTD., a Japanese corporation, ("NIDEK") for software (including but not limited to software linked dynamically or statically with other software) supplied by NIDEK or its designee pursuant to this Agreement, whether software alone or embedded software in a NIDEK hardware product, whether on disk or in read only memory, or on other media, or through an authorized website or network, and any accompanying documents or materials (including, but not limited to, operation manuals and electronic documents for such software, and other software for displaying or saving the data acquired from or through other NIDEK hardware product) (collectively, the "Software").

The Software and NIDEK hardware product (collectively, "NIDEK product") may include a third party's software which is linked, whether dynamically or statically, with the Software (the "Third-Party-Software"). The Third-Party-Software shall not be included in the definition of the "Software" in this Agreement. The rights and title of the Third-Party-Software belong to the third party, and the terms of use of the Third-Party-Software are set forth separately from this Agreement. The terms in this Agreement will not apply to the use of the Third-Party-Software except as expressly stipulated herein.

By using or installing the Software, you agree to be bound to the terms and conditions of this Agreement. If you do not agree with this Agreement, please do not use or install the Software and return the Software to the company from which you obtained the Software.

1. GRANT OF LICENSE

- 1.1. Subject to the terms and conditions set forth in this Agreement, NIDEK grants to you, and you accept, a limited, non-transferable and non-exclusive license to use the Software.
- 1.2. Unless otherwise agreed in writing by NIDEK or its designee, the license is limited to using the Software on a single computer or a single NIDEK hardware product and if you replace such computer or NIDEK hardware product, you may not use the Software without a new license of the Software.
- 1.3. Notwithstanding the provision of 1.2, if you connect a single server computer with the Software installed to a plurality of client computers, you may use the Software on such client computers; provided, however, that the upper limit of the number of said client computers will be determined by NIDEK in writing separately and individually from this Agreement.

- 1.4. Notwithstanding the provision of 1.2, if NIDEK permits you to install the Software on a plurality of computers using one license key of the Software, you may install and use the Software on such computers up to the upper limit of the number determined by NIDEK in writing separately and individually from this Agreement.
- 1.5. The Software is only to be used for its intended purpose provided in the specifications, operation manual or related documents in accordance with applicable laws and regulations. If the Software is embedded software in a NIDEK hardware product, you will use such Software only as embedded software for the use of such NIDEK hardware product.
- 1.6. For the license of the Software granted in this Agreement, unless the license is granted by NIDEK or its designee explicitly free of charge, you will pay to NIDEK or its designee the price for the Software, or if the Software is embedded software in a NIDEK hardware product, the price for the NIDEK hardware product in which the Software is embedded.

2. INTELLECTUAL PROPERTY RIGHTS

- 2.1. NIDEK, or an owner of the Third-Party-Software, retains any and all legal rights, title and interests in and to the Software or the Third-Party-Software. Any and all rights under copyright law, patent law, design law and other intellectual property laws not expressly granted herein are reserved by NIDEK or the owner of the Third-Party-Software. The license granted herein will not be intended as, or construed to be, any assignment of the rights of NIDEK or the owner of the Third-Party-Software. The Software and the Third-Party-Software are protected by copyright and other intellectual property laws and international treaties.

3. LIMITATIONS

- 3.1. You may not use the Software for any products without a license of the Software.
- 3.2. Unless otherwise permitted and other than the part specified by NIDEK in operation manuals or any accompanying documents for the Software, you may not analyze, reverse-engineer, decompile, disassemble or otherwise attempt to discover the source code of the Software.
- 3.3. You may not alter, reproduce, modify, translate, adapt, or divert the Software.
- 3.4. You may not remove, delete or change the copyright notice or other legends of the Software.
- 3.5. You may not sell, distribute, rent, license, sublicense, lease, assign or otherwise transfer the Software to third parties, or operate the Software for the benefit of third parties without prior written consent of NIDEK.
- 3.6. You may not create derivative works or cause or permit others to create derivative works based upon the Software without prior written consent of NIDEK.
- 3.7. You may not disclose operation manuals for the Software to any third party without prior written consent of NIDEK; provided, however, for the avoidance of doubt, the "third party" in this section will not include doctors, examiners, nurses, employees, patients and other persons who need to know the Software.
- 3.8. You may not use NIDEK's trademarks or trade names without prior written consent of NIDEK.

4. CONDITIONS OF USE

- 4.1. You shall take necessary measures (including but not limited to antivirus software) to prevent failure of NIDEK product due to external factors; provided, however, that in the case where it is otherwise provided in the provisions of operation manuals for NIDEK product or other documents, you shall take such necessary measures to the extent not inconsistent with such provisions.
- 4.2. If you enter data into NIDEK product or obtain data by the use of NIDEK product, you shall obtain and save backup of such data.

5. EXPORT RESTRICTIONS

- 5.1. If you export or re-export, directly or indirectly, the Software, you must comply with applicable export laws and regulations of Japan and other countries, and obtain any licenses or approvals required by governmental authorities.

6. UPDATES

- 6.1. The Software and/or the Third-Party-Software may be, at NIDEK's own discretion, changed, updated or modified from time to time without any prior notice to you. If such changes, updates, and modifications are applied to the Software licensed to you under this Agreement, such changes, updates, and modifications will be deemed a constituent part of the Software, and the terms and conditions of this Agreement will apply to such changes, updates, and modifications.
- 6.2. NIDEK may, at its own discretion, make amendments to any provisions of this Agreement (the "Amendments"), if NIDEK deems that:
 - a) such Amendments are appropriate in terms of interests for customers of this Software; or
 - b) such Amendments are commercially reasonable and not contrary to the objective of this Agreement, even if such Amendments are disadvantageous to you.Prior to the amendments, NIDEK will notify you of the terms and the effective date of such Amendments on the website or by any other means.
- 6.3. If you use the Software after the effective date of such Amendments, you shall be deemed to have agreed to such Amendments.

7. TERMINATION

- 7.1. This Agreement is effective until terminated. If you breach any term or condition of this Agreement, NIDEK may, without giving any prior notice to you, terminate this Agreement with immediate effect. Upon termination of this Agreement due to the breach of this Agreement, NIDEK reserves all the rights to claim damages resulting from such breach.
- 7.2. If this Agreement is terminated in accordance with the provision of 7.1., you must immediately cease the use of the Software, and delete, destroy and erase all the Software. Any fees paid by you for the license of the Software will not be refund for any reasons.

8. NO WARRANTIES

- 8.1. NIDEK MAKES NO REPRESENTATIONS OR WARRANTIES OF ANY KIND, EXPRESS OR IMPLIED, CONCERNING THE SOFTWARE AND THE THIRD-

PARTY-SOFTWARE, INCLUDING, WITHOUT LIMITATION, WARRANTIES OF MERCHANTABILITY, FITNESS FOR A PARTICULAR PURPOSE, NON-INFRINGEMENT OF THIRD PARTY RIGHTS, INCLUDING, WITHOUT LIMITATION, THIRD PARTY INTELLECTUAL PROPERTY RIGHTS, ACCURACY, RELIABILITY OR AVAILABILITY, ABSENCE OF OR RECOVERY FROM ANY INTERRUPTION, ERROR-FREE OPERATION OR CORRECTION OF DEFECTS OR MALFUNCTIONS.

9. LIMITATION OF LIABILITY

- 9.1. EXCEPT OTHERWISE EXPRESSLY STIPULATED IN THIS AGREEMENT, IN NO EVENT WILL NIDEK BE LIABLE FOR ANY INCIDENTAL, INDIRECT, SPECIAL, PUNITIVE, OR CONSEQUENTIAL DAMAGES, LOSS, CLAIMS OR COSTS WHATSOEVER, INCLUDING, WITHOUT LIMITATION, ANY LOST DATA, PROFITS, REVENUES, BUSINESS OPPORTUNITIES OR INFORMATION, LOSS OF USE OF ANY PRODUCT, PROPERTY OR EQUIPMENT, DOWNTIME COST, COST OF PROCUREMENT OF SUBSTITUTE GOODS OR SERVICES, OR ANY CLAIMS BY A THIRD PARTY, ARISING OUT OF OR RELATED TO THE USE OR INABILITY TO USE THE SOFTWARE AND/OR THE THIRD-PARTY-SOFTWARE, CHANGES, UPDATES OR MODIFICATIONS OF THE SOFTWARE AND/OR THE THIRD-PARTY-SOFTWARE, OR MAINTENANCE OR REPAIR SERVICE OF THE SOFTWARE IF ANY (collectively, the "DAMAGES"). THE ABOVE LIMITATIONS WILL APPLY REGARDLESS OF THE FORM OF ACTION, WHETHER IN CONTRACT, TORT, STRICT PRODUCT LIABILITY, OR OTHERWISE, EVEN IF NIDEK IS NOTIFIED OF THE POSSIBILITY OF SUCH DAMAGES.
- 9.2. THE LIMITATIONS PROVIDED IN THE PROVISION OF 9.1. SHALL NOT APPLY IN THE CASE WHERE THE DAMAGES ARE ATTRIBUTABLE TO NIDEK OR NIDEK IS LIABLE FOR SUCH DAMAGES IN ACCORDANCE WITH THE LAWS. EVEN IN SUCH CASE, NIDEK SHALL NOT BE LIABLE FOR ANY CONSEQUENTIAL, INDIRECT, INCIDENTAL, PUNITIVE OR SPECIAL LOSS OR DAMAGE. NIDEK'S TOTAL AGGREGATE LIABILITY FOR THE DAMAGES SHALL NOT EXCEED AN AMOUNT ACTUALLY PAID BY YOU FOR PURCHASE OF NIDEK PRODUCT; PROVIDED, HOWEVER, THAT THE LIMITATION OF THE AMOUNT SHALL NOT APPLY IN THE CASE WHERE THE APPLICABLE LAW PROHIBITS SUCH LIMITATION OR THE DAMAGES ARISING FROM NIDEK'S GROSS NEGLIGENCE OR WILLFUL MISCONDUCT.

10. GOVERNING LAW AND ARBITRATION

- 10.1. This Agreement will be governed by and construed in accordance with the laws of Japan.
- 10.2. All disputes arising between you and NIDEK relating to this Agreement or the interpretation or performance thereof will be finally settled by binding arbitration in Tokyo in accordance with the Commercial Arbitration Rules of The Japan Commercial Arbitration Association. Judgment upon the award rendered by arbitration will be final and may be entered in any court having jurisdiction thereof.

11. SEVERABILITY

11.1.If any provision or any portion of any provision of this Agreement will be held to be invalid or unenforceable, that provision will be severed from this Agreement and such invalidity or unenforceability will not affect the remaining provisions of this Agreement. The remaining provisions of this Agreement will continue in full force and effect.

12. SURVIVAL

12.1.The provisions of 2, 3, 5, 7, 8, 9, 10, 11, 13, 14, 15, 16, 17, 18, 19 and this provision will survive the termination of this Agreement and will be binding after the termination of the Agreement.

13. ASSIGNMENT

13.1.This Agreement or any part of this Agreement may not be assigned or transferred without prior written consent of NIDEK. The permitted assignee or transferee must agree to all the terms and conditions of this Agreement prior to the assignment or transfer.

13.2.This Agreement will be binding upon the permitted assignee or transferee and be enforceable by NIDEK.

14. ENTIRE AGREEMENT

14.1.This Agreement constitutes the entire agreement between you and NIDEK concerning the Software, and supersedes any prior written or oral agreement between you and NIDEK. No modification of this Agreement will be binding unless otherwise agreed in writing.

15. NO WAIVER

15.1.The failure of NIDEK to enforce at any time or for any period the provisions hereof in accordance with its terms will not be construed to be a waiver of such provisions or of the rights thereafter to enforce each and every provision.

16. NO THIRD PARTY RIGHTS

16.1.This Agreement is intended to be solely for the benefit of you and NIDEK and is not intended to confer any benefits upon or create any rights in favor of any person other than you and NIDEK.

17. HEADINGS

17.1.All headings are for convenience only and will not affect the meaning of any provision of this Agreement.

18. LANGUAGE

18.1.The license agreement for the Software may be provided in multiple languages. In such event, unless otherwise agreed in writing, the following shall apply:

- a) If you use the Software in any countries outside Japan, the license agreement for the Software shall be executed and delivered in a text using the English language. The text using the English language shall prevail and control; and
- b) If you use the Software in Japan, the license agreement for the Software shall be executed and delivered in a text using Japanese language. The text using the Japanese language shall prevail and control.

19. APPLICATION OF SOFTWARE LICENSE AGREEMENT

19.1.If the terms and conditions of the "Software License Agreement" included in operations manuals for NIDEK product are inconsistent with the terms and conditions of the "Software License Agreement" displayed on NIDEK product, the terms and conditions of the "Software License Agreement" included in operations manuals for NIDEK product prevail.

Table of Contents

1 SAFETY PRECAUTIONS - - - 7

- 1.1 For Safe Use - - - 7
- 1.2 Signal Words for Safety - - - 7
- 1.3 Cautions for Administration and Controlled Area - - - 8
- 1.4 Usage Precautions - - - 9
- 1.5 Labels and Symbols - - - 14

2 INTRODUCTION - - - 17

- 2.1 Outline - - - 17
- 2.2 Intended Use - - - 17
- 2.3 Intended Patient Population - - - 17
- 2.4 Intended User Profile - - - 17
- 2.5 Intended Use Environment - - - 18
- 2.6 Principles - - - 19
- 2.7 Precautions in Patient Selection - - - 20
- 2.8 Adverse Events and Adverse Device Effects - - - 20
- 2.9 Packed Contents - - - 21
- 2.10 Device Configuration - - - 23
 - 2.10.1 Connector box - - - 26
 - 2.10.2 Control box - - - 27
 - 2.10.3 Optional optical table - - - 29
- 2.11 Screens and Functions - - - 30
 - 2.11.1 Main screen (YAG mode) - - - 30
 - 2.11.2 Main screen (SLT mode) - - - 33
 - 2.11.3 Summary of treatment (YAG mode) - - - 36
 - 2.11.4 Summary of treatment (SLT mode) - - - 37

3 OPERATING PROCEDURE - - - 39

- 3.1 Operation Flow - - - 39
- 3.2 Operating Procedure - - - 41
 - 3.2.1 Activating the device - - - 41
 - 3.2.2 Preparing for emission - - - 42
 - 3.2.3 Emitting the treatment beam - - - 45
 - 3.2.4 Stopping the device - - - 47
- 3.3 Using SLT-NAVI (SLT Mode) - - - 48
- 3.4 Changing Device Settings - - - 50
 - 3.4.1 Changing the settings - - - 50
 - 3.4.2 Setting items - - - 51
- 3.5 Using Limit Mode (YAG Mode) - - - 55
 - 3.5.1 Using Limit mode from start-up - - - 55
 - 3.5.2 Using Limit mode after activating YAG mode - - - 56
 - 3.5.3 Deactivating Limit mode - - - 58

-
- 3.6 Using Memory - - - 59
 - 3.6.1 Loading memory - - - 59
 - 3.6.2 Registering memory - - - 60
 - 3.6.3 Deleting memory - - - 62

4 MAINTENANCE - - - 65

- 4.1 Troubleshooting - - - 65
- 4.2 Error Messages and Remedies - - - 67
- 4.3 Replacing Consumables - - - 69
 - 4.3.1 Consumable list - - - 69
 - 4.3.2 Attaching chinrest paper - - - 69
 - 4.3.3 Replacing fuses of optional optical table - - - 70
- 4.4 Cleaning - - - 71
 - 4.4.1 Cleaning the device exterior - - - 71
 - 4.4.2 Cleaning areas that come into contact with patients - - - 71
 - 4.4.3 Cleaning optical parts - - - 72
- 4.5 Checking Optical Axis of Treatment Beam (YAG mode) - - - 73
- 4.6 Checking Optical Axis of Treatment Beam (SLT mode) - - - 75
- 4.7 Measuring and Calibrating Treatment Beam Energy Output - - - 77
 - 4.7.1 Measuring treatment beam energy output - - - 77
 - 4.7.2 Calibrating treatment beam energy output - - - 79
- 4.8 Measuring and Calibrating Aiming Beam Power Output - - - 84
 - 4.8.1 Measuring aiming beam power output - - - 84
 - 4.8.2 Calibrating aiming beam power output - - - 86

5 SPECIFICATIONS AND TECHNICAL INFORMATION - - - 87

- 5.1 Specifications - - - 87
- 5.2 Light Hazard - - - 92
- 5.3 EMC (Electromagnetic Compatibility) - - - 94
- 5.4 Administration List - - - 97
- 5.5 Laser Beam Output Calibration Record - - - 98
- 5.6 Pre-use Checklist - - - 100
- 5.7 Glossary - - - 101

6 INDEX - - - 105



SAFETY PRECAUTIONS



1.1 For Safe Use

 BEFORE USE, READ THIS MANUAL.
Be sure to read the operator's manual prior to use of the device to understand the safety precautions and operating procedures thoroughly. Keep this manual handy for reference.

1

1.2 Signal Words for Safety

In this manual, signal words are used to designate the degree or level of safety alerting. The definitions are as follows.
 WARNING Indicates a potentially hazardous situation which, if not avoided, may result in death or serious injury.
 CAUTION Indicates a potentially hazardous situation which, if not avoided, may result in minor or moderate injury or property damage accident.
Even situations indicated by  CAUTION may result in serious injury under certain conditions. Safety precautions must be strictly followed at all times.

1.3 Cautions for Administration and Controlled Area

Administration

CAUTION

- The medical institution in which the device is to be used must appoint at least one qualified administrator and one sub-administrator.
- Administrators are responsible for the storage and administration of the device within the area in which the device is to be used.
- Administrators should appoint personnel who can use the device, conduct necessary training, and provide them with the latest information (seminars, workshops, and conferences may be regarded as a training).
- Authorized personnel must follow the instructions of the administrators.
- Administrators should prepare and keep a registered user list.
[↪ “5.4 Administration List” \(page 97\)](#)
- Authorized personnel must be appointed by the administrators and are required to thoroughly understand the operating procedures and safety control of the device.

Controlled area

CAUTION

- A controlled area in which the device is to be used must be specified by the medical institution. Indications specifying the controlled area are required (Controlled area indication).
Attach the DANGER label to the entrance of the controlled area.
[↪ “2.9 Packed Contents” \(page 21\)](#)
- In the controlled area, notices that are necessary for administration of the device such as the laser name, warnings, and other information should be indicated (Warning indications).
- Any persons entering the controlled area (except the authorized personnel specified in the registered user list) are required to obtain permission from the administrators, understand the precautions, and take necessary protections prior to entering the controlled area.
- It is recommended that all personnel entering the controlled area take a visual acuity test before entering and after leaving the controlled area to make sure that their visual acuity is not impaired.

Preparation of controlled area facilities and equipment

CAUTION

- Administrators should prepare the facilities necessary for introducing the device.
- Administrators should prepare the facilities and all equipment necessary for maintenance and safety control of the device.
- Administrators should regularly conduct maintenance and checks as described in this operator’s manual, and record the results.

1.4 Usage Precautions

Installation

WARNING

- Only service personnel trained by NIDEK can install and adjust the device.
- Ensure that there is no flammable anesthetic gas in the operating room.
Laser emission may cause fire or explosion.
- Be sure to connect the power plug to a grounded power outlet.
Electric shock or fire may result from device malfunction or electric leakage.

CAUTION

- Install the device in an environment that meets the specifications.
 *"5.1 Specifications" (page 87)*
- Install the device in a location that is not exposed to water or rain.
Malfunction may result.
- Be sure to use a (hospital grade) power outlet which meets the power specification requirements.
The device may not work properly, or malfunction or fire may result.
- Do not use a power cord other than the one provided. Also do not connect the provided power cord to any other device.
- Install the device so that the outlet that the power plug is inserted into is easily accessible during use. In addition, ensure that the power plug can be disconnected without the use of any tool.
Failure to do so may interfere with disconnecting the power from the input power source in case of an abnormality.
- Securely connect the power cord and cables to the specified connectors.
Malfunction may result.
- Do not coil the power cord or cables too tightly, or crush or pinch them with heavy objects.
Electric shock or fire may result.
- If the internal wires are exposed, replace the cords and cables with new ones.
Electric shock or fire may result.
- Immediately replace the cords and cables if the power is intermittent, or the cord or plug is hot to touch.
Malfunction or fire may result.
- Do not drag the device by the power cord or cables.
Cables may break and the device may topple causing injury or malfunction.

Before use

WARNING

- If any serious device-related incident occurs, report it to NIDEK and the competent authority in the country where the user or patient, or both reside.
- Do not use this device for any purposes other than the intended purpose.
 *"2.2 Intended Use" (page 17)*
- Use of this device is limited to ophthalmologists.
Unanticipated adverse events and adverse device effects may result.
- Be sure to read the operator's manual prior to operation of the device to understand the safety precautions and operating procedures thoroughly.
Use of the device outside the scope may cause unanticipated adverse events and adverse device effects.
- Do not perform operations by procedures different from those described in the operator's manual.
Control and adjustment of the device in procedures outside the scope may cause exposure to hazardous laser radiation.
- Use the device with at least one assistant in the room.
This is in precaution to a case such as electric shock. It is desirable that the assistant is trained in resuscitation.
- Backup measures for the scheduled surgery must be prepared.
This is in precaution to device failure.
- Pay attention when using the device with other equipment that comes into contact with the patient.
Electromagnetic wave or other interference may cause a danger.
Using an electrocautery for coagulation may cause electric shock or burns to the contact area.
- Do not use accessories other than those specified by NIDEK.
Use of the device outside the scope may cause unanticipated adverse events and adverse device effects.
- Do not modify or touch the device interior.
Electric shock or malfunction may result.
- Have personnel in the operating room wear safety goggles. Instruct them never to gaze directly at the treatment beam even while wearing the safety goggles.
YAG treatment beam: Wavelength 1,064 nm, OD > 7, 950-1080 D LB6 + IRM LB7 (EN207)
SLT treatment beam: Wavelength 532 nm, OD > 7, 315-532 DIRM LB6 (EN207)
- When the device is not in use, remove and store the keycard in a secure place.
Accidents or information leakage caused by unauthorized personnel may result.
- When connecting to other devices via a networked medical system, confirm that no hazard affects patients, operators, or third parties. In addition, when a device is added or removed from the network, or the device is updated or upgraded, confirm that no hazard affects patients, operators, and third parties in the same manner.

During use

WARNING

- Before using the device, perform operation check and record the results. If abnormal conditions are encountered during the check, do not use the device.
If the device is used under abnormal conditions, intended results may not be achieved. Also unanticipated malfunctions or health hazards due to improper diagnosis may occur.
- Check the setting values before emitting the treatment beam.
The treatment beam may be emitted in unintended settings.
- Start the treatment beam from the lowest energy output, then gradually increase the energy until the desired effect is obtained. Be sure to return the output energy to the lowest after every operation.
Excessively intense treatment beam may be emitted.
- When laser beam emission is not intended (such as during observation of an eye), set the device to STANDBY mode so that laser emission is not possible.
Accidental exposure to the laser beam may result.
- Before connecting and removing the power cord from the power outlet, confirm that the master switch is turned off.
Electric shock or malfunction may result.
- If any abnormal indication (other than laser beam emission conditions) is displayed on the control box, follow the applicable instructions.
 [“4.2 Error Messages and Remedies” \(page 67\)](#)
- Before lowering the table, check that there are no objects that may get caught in it.
The feet of the patient or operator may be caught and may be injured.

CAUTION

- Do not perform servicing or maintenance on the device during use.
- Prior to surgery, provide the patient with sufficient information about the surgery purposes and methods, expected results, and possible adverse events.
- If the illumination light is not perceived, stop the surgery.
The affected area may not be observed. Unanticipated adverse effects on laser beam emission may result.
- Always pay attention to the direction of the aiming beam.
When the treatment beam is applied to tissue, damage to the cornea and such, blindness, pain of skin, or burn injury may result.
- Do not look directly at the aiming beam or point the beam toward others.
Unnecessary exposure to the laser beam may result.
- Before emitting the treatment beam, confirm that there is no reflective object in the laser optical path.
Exposure to the reflected laser beam may result.
- Take care not to catch hands or fingers in moving parts. Be sure to also caution patients.
Hands or fingers may be pinched and may result in injury.
- When the patient comes off from the device, instruct the patient not to stand up while holding the chinrest.
The device may topple over resulting in injury.
- If the device fails, disconnect the power cord from the power outlet without touching the device interior.
Contact NIDEK or your authorized distributor.


CAUTION

- If power failure or device failure occurs while the patient is wearing the head belt, loosen the head belt immediately.

Unanticipated malfunctions or adverse events may result.

- Use devices that comply with IEC 60601-1 in the patient environment. If any device that does not comply with IEC 60601-1 is to be used, use an isolation transformer or common protective grounding.

The patient environment is the volume of space in which contact can occur between the patient and any part of the device (including connecting devices) or between the patient and any other person(s) touching the device (including connecting devices).

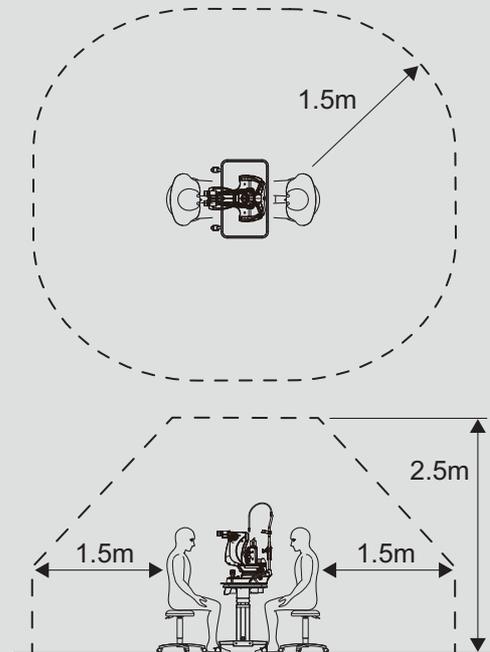

After use

CAUTION

- Clean the forehead rest, chinrest, grips, and head belt for each patient. Use alcoholic disinfectants if necessary.

When using the chinrest paper, remove one sheet of paper.

- After using the device, turn it off and cover it with its dust cover.

Dust may affect the performance of laser beam emission.

- After using the device, disconnect the power cord from the power outlet.

During transport

CAUTION

- Before moving the device to another room, tighten the knobs to secure the device.

Hands or fingers may be pinched and may result in injury.

- Store the device in an environment that meets the specifications during transport and storage.

 "5.1 Specifications" (page 87)

- Before transporting the device, pack it in the packaging material.

Pack it with the moving parts locked and the slit opened.

Malfunction may result.

- Do not bump the device during transport.

Optical axis may be shifted.

 **CAUTION**

- Care should be taken so that the temperature varies as little as possible during transport. Changes in temperature may result in condensation inside the device or adversely affect the functionalities.
-
-

Maintenance

 **CAUTION**


- Only service personnel trained by NIDEK can repair the device.
NIDEK assumes no responsibility for any adverse events resulting from improper servicing.
 - To ensure the continued safe use of the device, the manager of this device must make sure that maintenance and preventive inspection are performed at least once a year.
For details of maintenance and preventive inspection, ask NIDEK or your authorized distributor. If the manager of this device cannot perform the maintenance and preventive inspection, contact NIDEK or your authorized distributor.
 - When sending the device back to NIDEK for repair or maintenance, clean the surfaces of the device (especially, the areas that come into contact with the patient) with a clean cloth dampened with rubbing alcohol.
 - It is recommended to replace the device before its service life expires.
Even with proper maintenance and check, after the expected service life, the device reliability or safety may become degraded and fail to achieve the target values.
-
-

Disposal

 **CAUTION**

- When disposing of the device, contact NIDEK or your authorized distributor.
 - Observe the local ordinances and recycling plans concerning disposal and recycling. Especially when disposing of lithium battery, printed circuit board, plastic parts containing brominated flame retardant, LCD, and power cord used internally, observe the instructions of local governments.
For details on local ordinances, contact your local governments.
 - When disposing of packaging materials, sort them by material and dispose according to the local ordinances and recycling plans.
-
-

1.5 Labels and Symbols

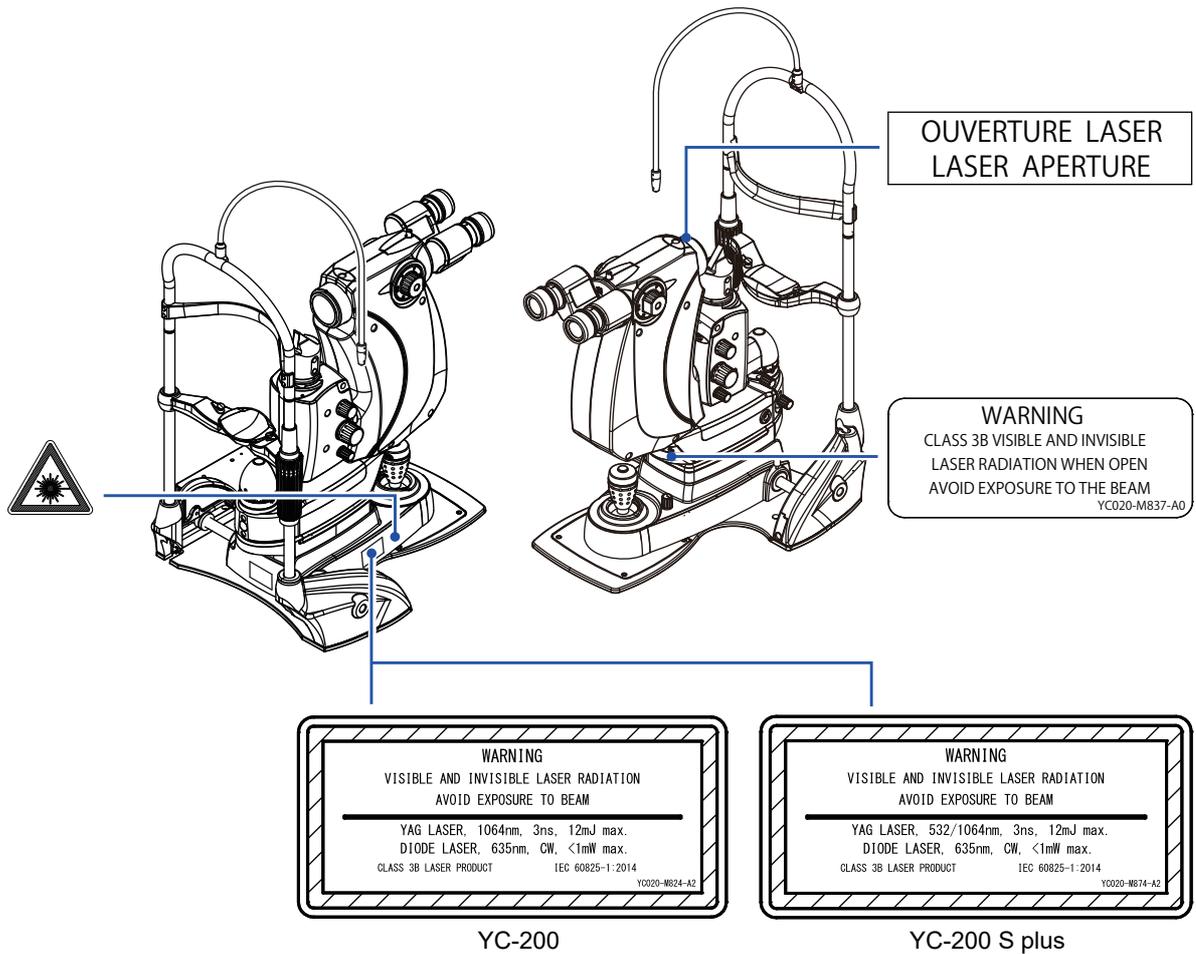

Refer to the operator's manual.



Type B applied part
The applied parts of the device are the forehead rest and chinrest.



“OFF” (power)



“ON” (power)



Alternating current



Start (of action)



Emergency laser stop



Warning; Laser beam
Background color: Yellow



Manufacturer



Date of manufacture



Indicates that this product is to be disposed of in a separate collection of electrical and electronic equipment in EU.



Medical device



EU authorized representative



Serial number



Unique device identifier



Catalog number



Swiss authorized representative

2

INTRODUCTION

2.1 Outline

This device is an ophthalmic pulsed laser system using a 1,064 nm Q-switched pulsed Nd: YAG laser as the treatment beam source.

Two types are available, differing in the available types of laser emission.*1

The operation mode available differs depending on the type.

Type name	Emitted laser (wavelength)	Operation mode available
YC-200	Nd: YAG laser (1,064 nm)	YAG mode(<i>page 30</i>)
YC-200 S plus	Nd: YAG laser (1,064 nm) SLT laser (532 nm)	YAG mode(<i>page 30</i>) SLT mode(<i>page 33</i>)

2.2 Intended Use

The OPHTHALMIC YAG LASER SYSTEM YC-200 is indicated for the performance of posterior capsulotomy, pupillary membranectomy, iridotomy (hole in the iris) and selective laser trabeculoplasty in phakic, aphakic and pseudophakic subjects.

2.3 Intended Patient Population

- Age
Infants to elderly (for whom the chinrest can be used)
- Health condition
Patient with eye disease requiring any of the treatments described in the intended use
- Conditions - Visual function
One or both eyes have disease.

2.4 Intended User Profile

Ophthalmologists

*1. This operator's manual describes the device using the YC-200 illustration.

2.5 Intended Use Environment

Medical facility

 **CAUTION**

- If the device is used outside the specified use location, intended performance and security level cannot be maintained.
-
-

2.6 Principles

This device uses passive Q-switching for laser oscillation. With passive Q-switching, the device first accumulates sufficient energy within the laser cavity then emits laser pulses of high peak power for extremely short durations by use of the Q-switched optical element which changes between transparent and opaque according to the accumulation of energy.

- YAG mode

The emitted YAG laser beam becomes coaxial with the aiming beam in the optical system within the laser system and converges at the intraocular target such as the posterior lens capsule after being emitted from the aperture of the device toward the patient's eye. When the YAG laser beam temporarily converges at a spot, and the power density exceeds a certain degree, plasma is formed. Sudden expansion and heating of the plasma generates shock waves. The shock waves destroy fine tissue around the focal point. In addition, the formed plasma absorbs and scatters the incident light. This feature protects the tissue posterior to the focal point.

A red diode laser (wavelength: 635 ± 5 nm) is used for the aiming beam that indicates the position to which the YAG laser beam is to be emitted. The device was designed so that the aiming beam becomes coaxial with the YAG laser beam in the optical system in the laser system, and that the spot where the two alignment spots coincide becomes the target of the YAG laser beam.

- SLT mode (YC-200 S plus)

The emitted YAG laser beam is converted to a 532 nm green laser beam (hereafter SLT laser beam) by means of a wavelength transducer within the device. In addition, it becomes coaxial with the aiming beam in the optical system within the laser system and is irradiated so that the spot diameter becomes 400 μm at the intraocular target such as open anterior chamber angle after being emitted from the aperture of the device toward the patient's eye. Selectively emitting the SLT laser beam onto the anterior chamber angle affects pigmented cells in the tissue to lower the intraocular pressure without destroying the affected tissue (less-invasive).

A red diode laser (wavelength: 635 ± 5 nm) is used for the aiming beam that indicates the position to which the SLT laser beam is to be emitted. The device was designed so that the aiming beam becomes coaxial with the SLT laser beam in the optical system in the laser system, and that the laser beam converges with the 400 μm laser emission size at the target.

2.7 Precautions in Patient Selection

Caution should be taken when using the device in regards to patients with the following conditions:

- Posterior capsulotomy when the IOL has become strongly adhered to the posterior capsule
- Patients with progressive eye disease
- Patients who have difficulty in eye fixation due to nystagmus or have a condition that may induce nystagmus
- Patients with opacity of the aqueous humor or cornea
- Patients whose intraocular lens is not sufficiently fixed after the surgery due to posterior capsule rupture or zonular rupture
- Patients with acute attack of primary angle closure (with corneal edema)
- Patients with advanced glaucoma with progressed visual field loss

2.8 Adverse Events and Adverse Device Effects

Possible adverse events and adverse device effects may include, but are not limited to the following:

● Adverse device effects

If any abnormality is found with the device during the pre-use check, do not use the device.

If the device fails, intended treatment results may not be obtained, and health hazards or unexpected adverse events described in [Adverse events] below may result.






● Adverse events





Possible adverse events (complications) may include, but are not limited to the following:

- Increased intraocular pressure
- Hyphema
- Bullous keratopathy
- Localized cataract
- Closure of perforating wound
- Pupillary block
- Iris injury
- Intraocular lens malpositioning and intraocular lens drop into the vitreous chamber
- Corneal burn
- Corneal endothelial damage
- Retinal detachment
- Ache
- Cystoid macular edema
- Anterior chamber inflammation
- Redness
- Corectopia
- Corneal opacity
- Irritation such as iritis, hyalitis, and uveitis
- Posterior synechia
- Peripheral anterior synechia (PAS)
- Corneal edema
- Damage to the intraocular lens
- Effects of erroneous irradiation to the retina
- Scleral perforation
- Macular hole
- Endophthalmitis
- Conjunctival hyperemia
- Anterior chamber angle hemorrhage
- Choroidal effusion
- Eye pain


2.9 Packed Contents

- The following are included in the standard configuration. Check the contents before use.

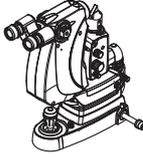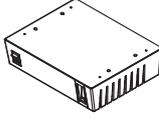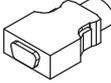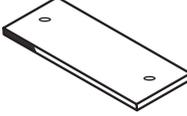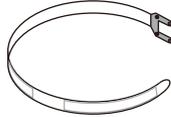	Part name	Quantity	Appearance
<input type="checkbox"/>	Main body	1 unit	
<input type="checkbox"/>	Head rest	1 unit	
<input type="checkbox"/>	Control box	1 unit	
<input type="checkbox"/>	Connector box	1 unit	
<input type="checkbox"/>	Remote connector plug	1 unit	
<input type="checkbox"/>	Keycard	1 unit	
<input type="checkbox"/>	Power cord	1 unit	
<input type="checkbox"/>	Chinrest paper	2 units	
<input type="checkbox"/>	Head belt	1 unit	

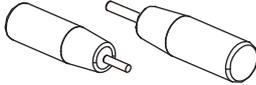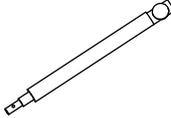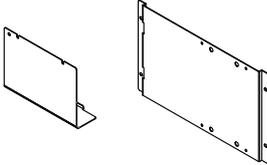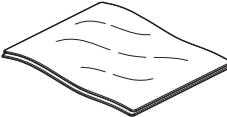	Part name	Quantity	Appearance
<input type="checkbox"/>	Grips	1 set	
<input type="checkbox"/>	Focusing rod	1 unit	
<input type="checkbox"/>	Attachment plates	1 set	
<input type="checkbox"/>	Arm rest	1 unit	
<input type="checkbox"/>	Cap	1 unit	
<input type="checkbox"/>	Dust cover	1 unit	
<input type="checkbox"/>	DANGER label	1 unit	
<input type="checkbox"/>	Operator's Manual	1 volume	

2.10 Device Configuration

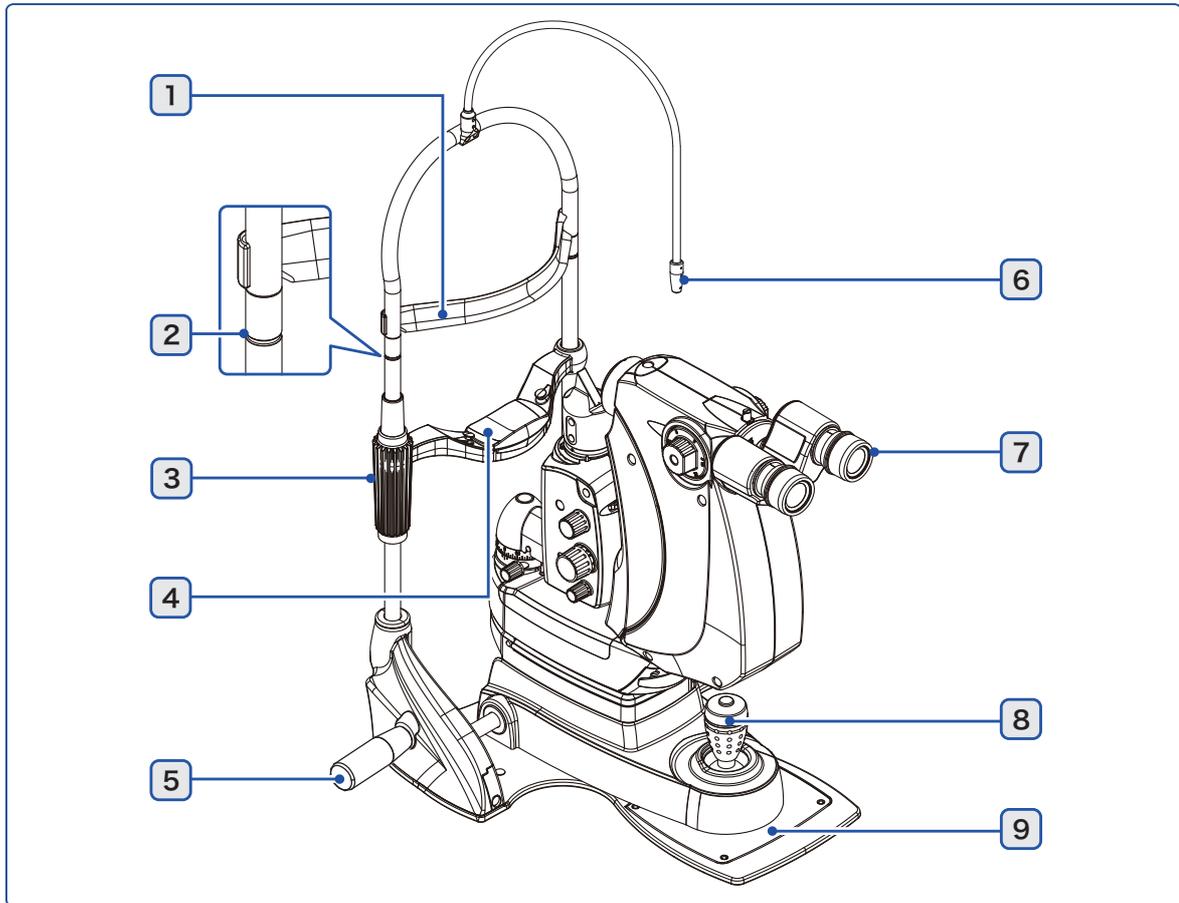

1 Forehead rest

2 Eye level marker

The patient's eyes are aligned to this height. Align the position by adjusting the height of the chinrest.

3 Chinrest elevation control

Adjusts the height of the patient's chin.

4 Chinrest

5 Grips

Patients hold them.

6 Fixation lamp

Used to steady the patient's visual axis by having the patient focus on the lamp. The position of the fixation lamp can be adjusted by the flexible arm.


7 Microscope

8 Joystick

Used for alignment and focusing.

By manipulating the joystick right and left, adjust right and left movement. By rotating the joystick, adjust up and down movement. Adjust the focus with forward and backward movements.

9 Slide plate

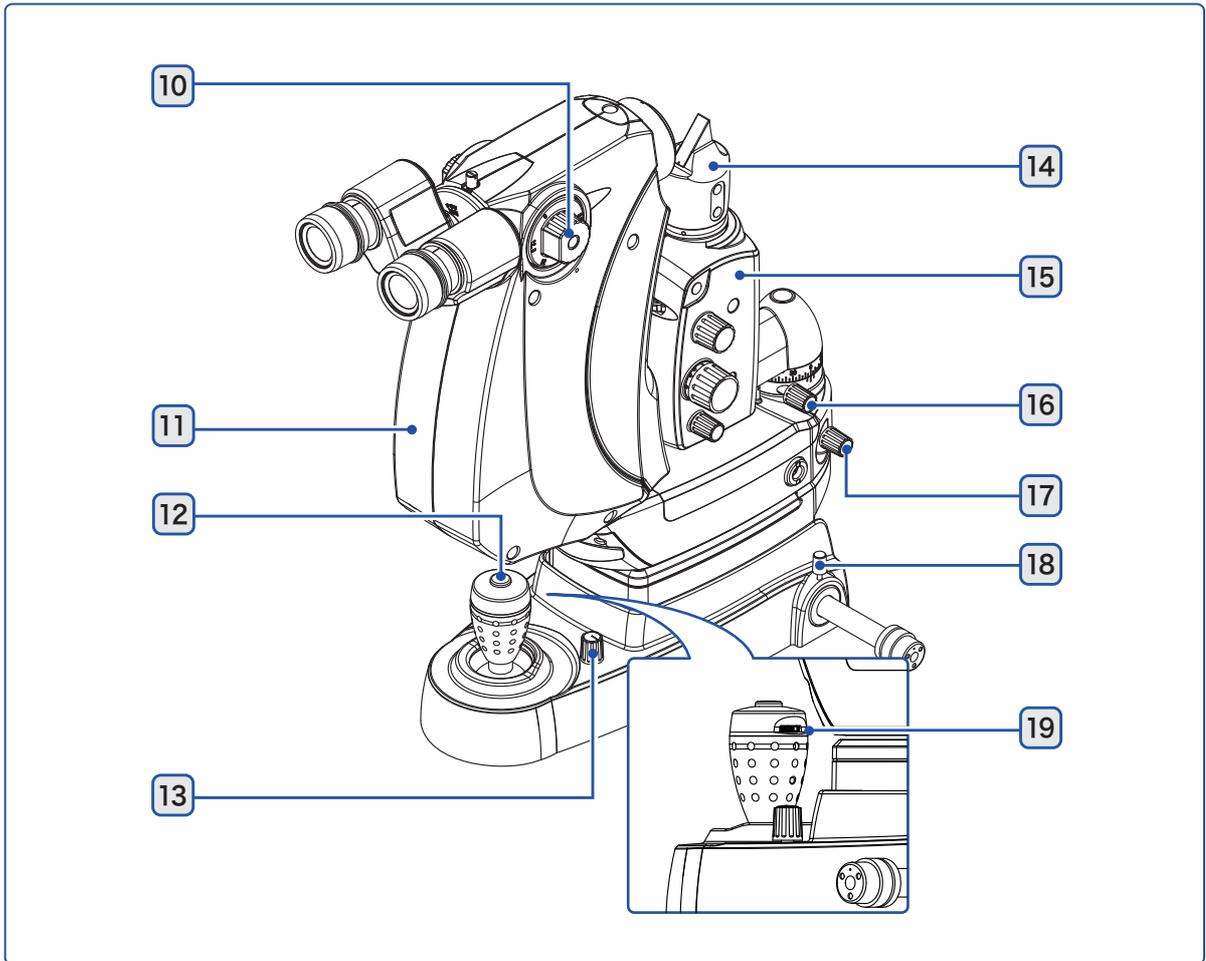

10 Magnification changer

Used to select the total magnification of the microscope.

11 Microscope arm

12 Hand switch


Used to emit the treatment beam. The beam is not emitted in STANDBY mode.

13 Illumination control

Used to adjust the intensity of the illumination light.


14 Illumination tower

Emits the illumination light rays that became parallel in the illumination optical system onto the affected area.


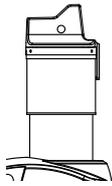
	
<p>YC-200</p>	<p>YC-200 S plus Equipped with a function to change the upward angle of the illumination light (tilting function) (Rotating the illumination tower changes the upward angle.)</p>

15 Illumination unit arm**16 Illumination unit arm fastening knob****17 Microscope arm fastening knob****18 Base unit fastening knob**

Used to fasten the device so that the device may not move horizontally.

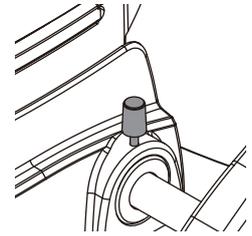**19 S-Switch**

An auxiliary switch built in the joystick. Three functions can be assigned as desired.

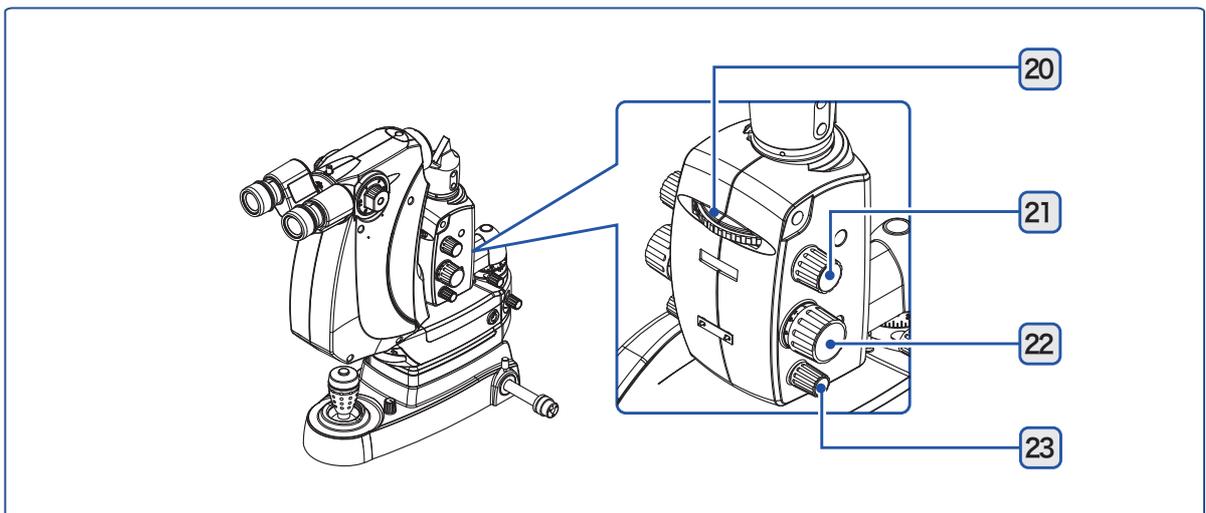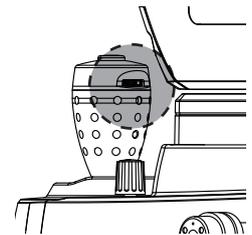**20 Filter changer**

Used to insert a filter into the illumination optical path.

21 Slit rotation control

Used to rotate the slit. A click is felt in the horizontal and vertical axes.


22 Slit width control

Used to adjust the slit width of the illumination light.

23 Slit length control

Used to adjust the aperture diameter or slit length of the illumination light.

2.10.1 Connector box

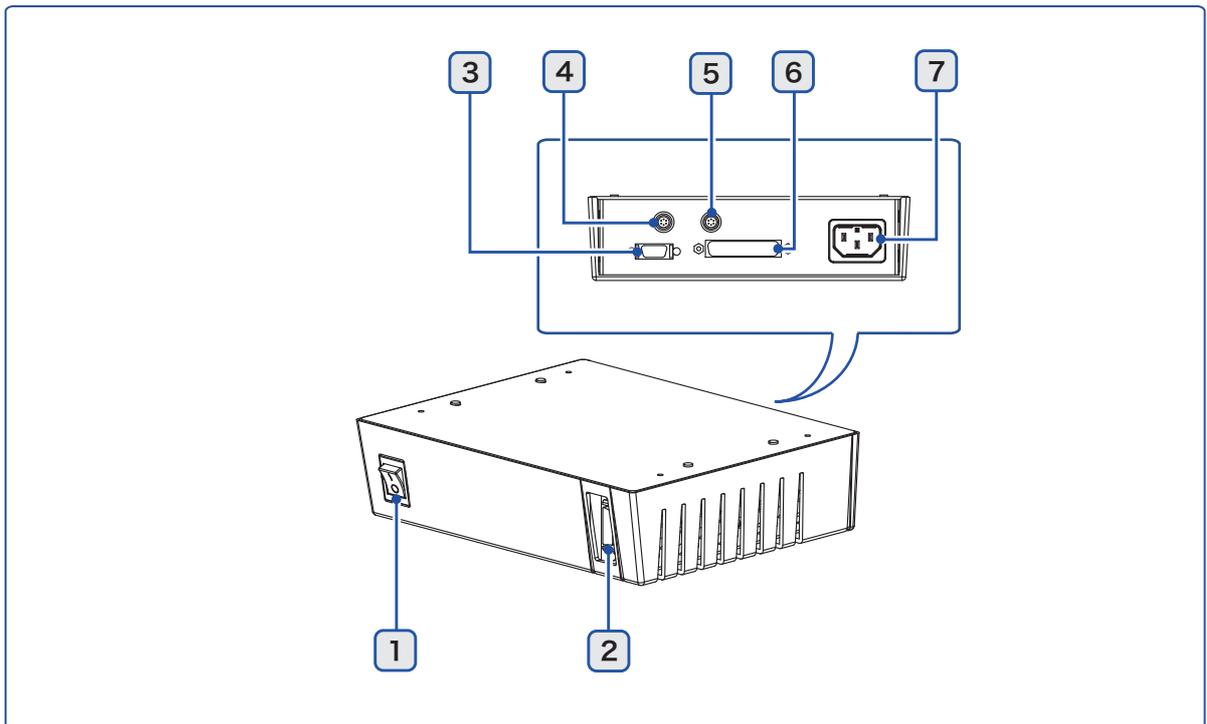

1 Master switch

2 CB connector

The control box cable is connected here.

3 REMOTE connector

An external switch for the remote interlock is connected here. If no external switch is to be used, insert the remote connector plug.

4 Foot switch connector

The optional foot switch cable is connected here.

5 Fixation lamp connector

The fixation lamp cable is connected here.


6 Main body connector

The main body cable is connected here.

7 Power inlet

The power cord is connected here.

2.10.2 Control box

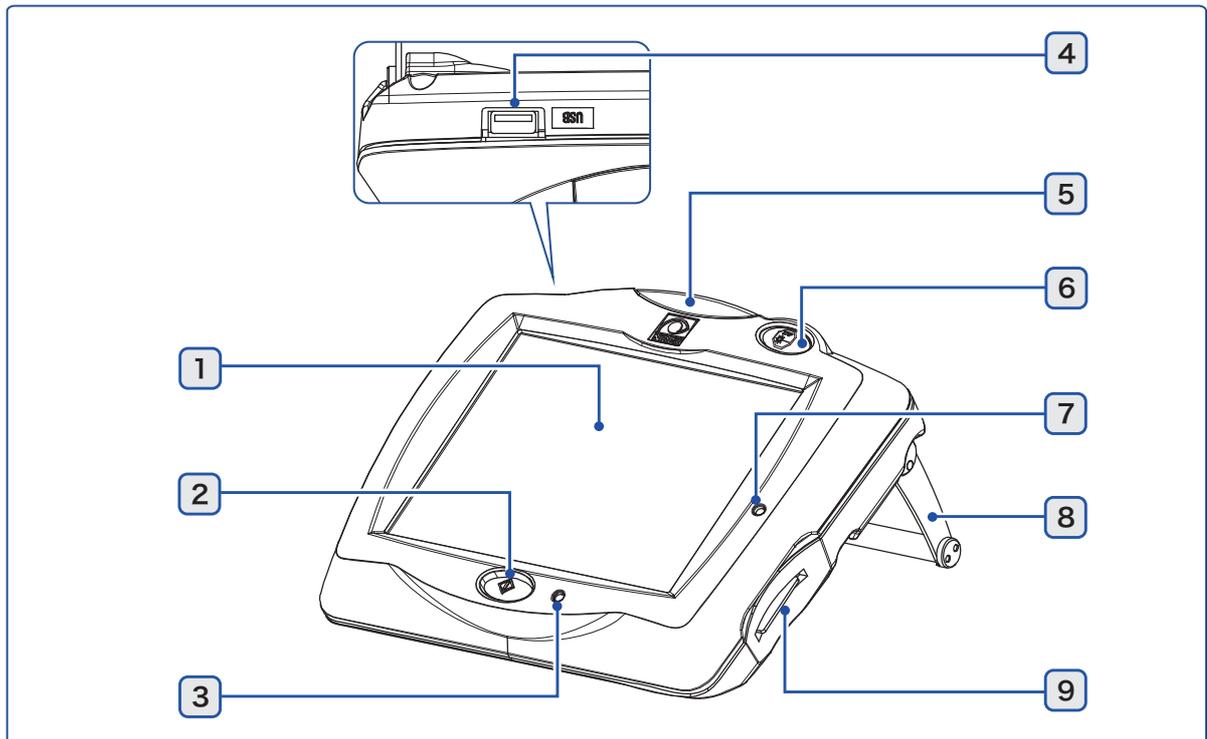

1 Touch screen

2 Start button

Used to activate or stop the device.

3 Power indicator

Lights up when the master switch is turned on (|).

4 USB-A port

Cannot be used.

5 Status indicator

Indicates the device status.

Status indicator	Status
Lit in orange	STANDBY (Treatment beam cannot be emitted.)
Lit in green	READY (Treatment beam can be emitted.)
Slow orange blink	Sleep mode ^{*a}

^{*a}. The device automatically enters sleep mode (the screen turns off) when the device is left idle for 10 minutes without any key operation. To return from sleep mode, touch the touch screen or press the start button of the control box.

6 Emergency laser stop button

Used to stop the device in case of an emergency.

To restart the device, press the start button.

7 Access indicator


Lights up when the keycard is inserted. The indicator blinks while data is being written to the keycard.

8 Stand


9 Keycard slot

The keycard is inserted here.

The device does not start without the keycard inserted.

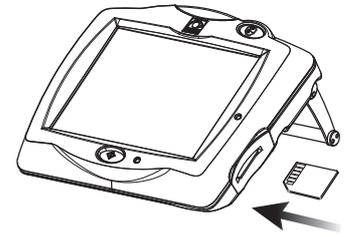

2.10.3 Optional optical table

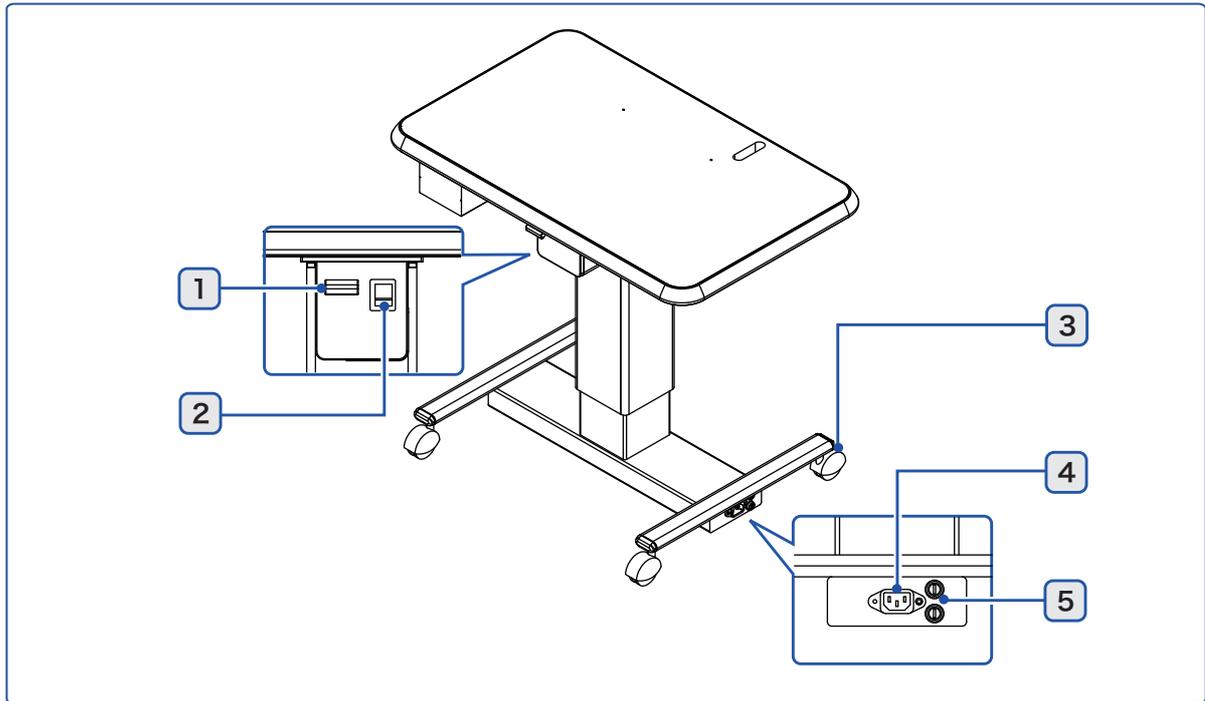

1 Table up/down switch

Used to move the optional optical table up and down.

2 Master switch

Turns the power on or off.

3 Caster lock

CAUTION

- Unlock the caster before moving the optional optical table.
The device may topple over resulting in injury.

4 Power inlet


The power cord is connected here.

5 Fuse holder^{*1}

^{*1}. The position of the fuse holder may differ depending on the destination.

2.11 Screens and Functions

2.11.1 Main screen (YAG mode)

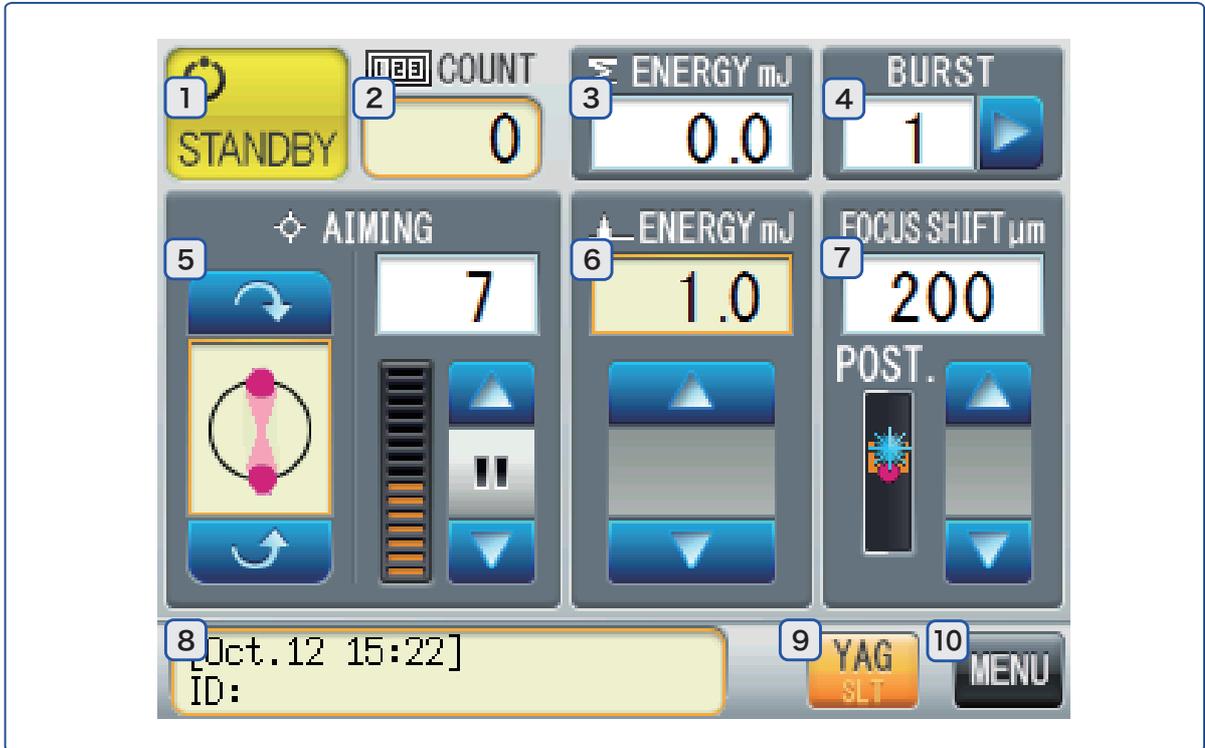

1 Status button

Toggles between READY (treatment beam can be emitted) and STANDBY (treatment beam cannot be emitted) modes.

Mode	Function
 (green)	To emit the treatment beam, set the device to READY mode. While the treatment beam is emitted, the EMISSION indicator  lights up in yellow.
 (yellow)	When the treatment beam is not to be emitted, set the device to STANDBY mode.


2 COUNT

Indicates the total number of treatment beam shots.



Knowledge

- Pressing the indication displays a confirmation window. Press "Yes" to reset the counter.

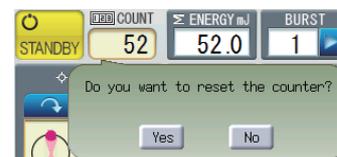

3 Σ ENERGY(TOTAL ENERGY)

Displays the integrated energy value of the emitted treatment beam.



Knowledge

- When the integrated energy is over 999.9 mJ, "OVER" is displayed.
- The value is cleared when the device power is turned on or the total number of treatment beam shots is cleared.
- The displayed integrated energy value is calculated from the setting values, not from the actual energy values.

4 BURST

Sets the number of pulses emitted in each laser emission.



5 AIMING

Sets the inclination and intensity of the aiming beam.

In conjunction with the set inclination, the inclination of  changes.




Button	Function
	Rotates clockwise.
	Rotates counterclockwise.
	Increases the intensity.
	Decreases the intensity.
	Turns on or off the aiming beam.



Knowledge

- Pressing  displays the pop-up window to select the angle from the specified ones.

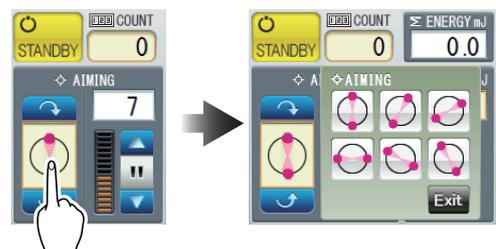

6 ENERGY


Sets the energy value of the treatment beam per pulse.

Button	Function
	Increases the energy output.
	Decreases the energy output.



Knowledge

- Pressing the indication displays the pop-up window to select the energy value from the predetermined options.

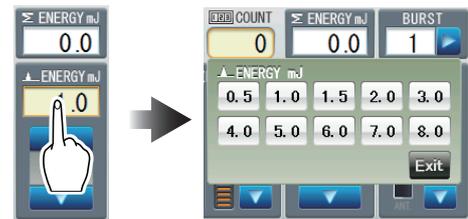

7 FOCUS SHIFT

Sets the focus position of the treatment beam relative to the aiming beam. The bar graph changes according to the specified focus position.

Button	Function
	Moves to the POST. (posterior chamber) side.
	Moves to the ANT. (anterior chamber) side.



8 Information display field

Displays various information. Contents can be selected by the parameter setting.

 "3.4 Changing Device Settings" (page 50)

9 Mode change button (YC-200 S plus)


Toggles between YAG mode and SLT mode.

10 MENU button

Various settings can be changed.

 "3.4 Changing Device Settings" (page 50)

2.11.2 Main screen (SLT mode)

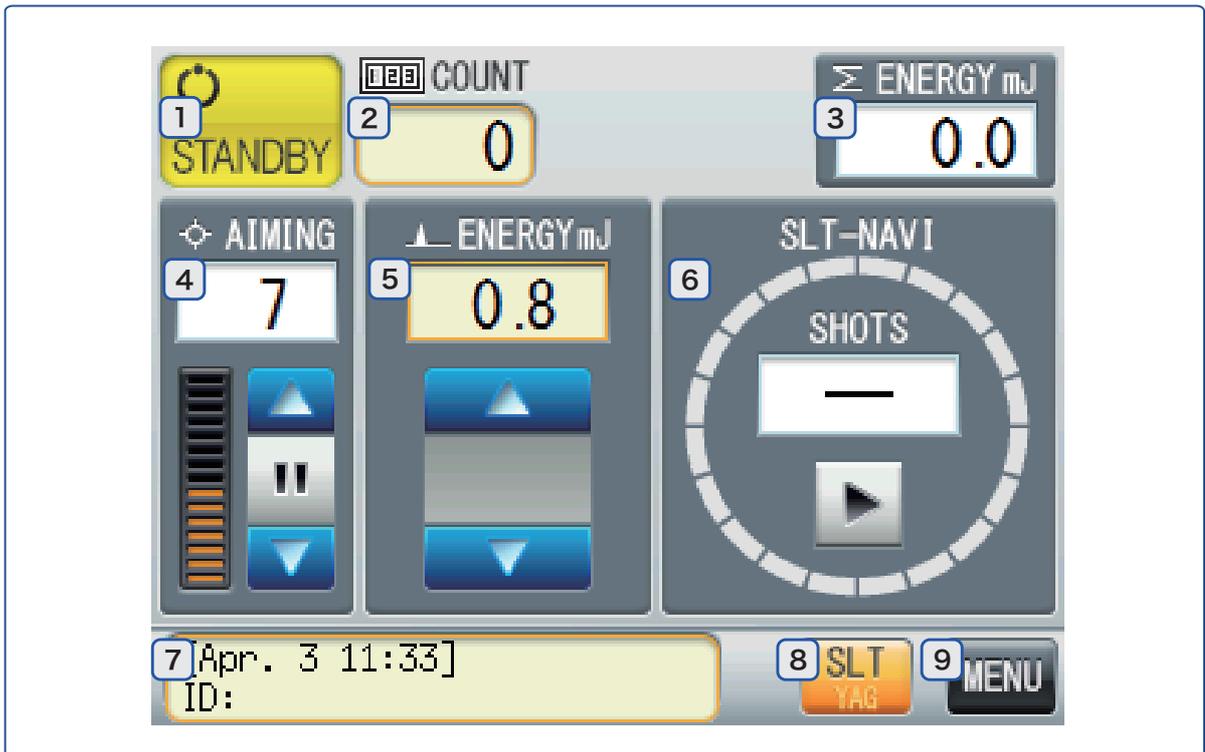

2

1 Status button

Toggles between READY (treatment beam can be emitted) and STANDBY (treatment beam cannot be emitted) modes.

Mode	Function
 (green)	To emit the treatment beam, set the device to READY mode. While the treatment beam is emitted, the EMISSION indicator  lights up in yellow.
 (yellow)	When the treatment beam is not to be emitted, set the device to STANDBY mode.

2 COUNT

Indicates the total number of treatment beam shots.



Knowledge

- Pressing the indication displays a confirmation window. Press "Yes" to reset the counter.



3 Σ ENERGY (TOTAL ENERGY)

Displays the integrated energy value of the emitted treatment beam.



Knowledge

- When the integrated energy is over 999.9 mJ, “OVER” is displayed.
- The value is cleared when the device power is turned on or the total number of treatment beam shots is cleared.
- The displayed integrated energy value is calculated from the setting values, not from the actual energy values.

4 AIMING

Sets the intensity of the aiming beam.

Button	Function
	Increases the intensity.
	Decreases the intensity.
	Turns on or off the aiming beam.



5 ENERGY


Sets the energy value of the treatment beam per pulse.

Button	Function
	Increases the energy output.
	Decreases the energy output.



Knowledge

- Pressing the indication displays the pop-up window to select the energy value from the predetermined options.

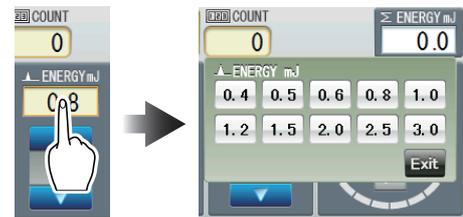

6 SLT-NAVI

Sets the emission area, emission start position, and emission direction.

↳ “3.3 Using SLT-NAVI (SLT Mode)” (page 48)

7 Information display field

Displays various information. Contents can be selected by the parameter setting.

↳ “3.4 Changing Device Settings” (page 50)



8 Mode change button 

Toggles between YAG mode and SLT mode.

9 MENU button 


Various settings can be changed.

 [“3.4 Changing Device Settings” \(page 50\)](#)

2.11.3 Summary of treatment (YAG mode)

The Summary of treatment screen is displayed at the time specified in “Summary Display Settings”.

↳ “3.4.2 Setting items” (page 51)

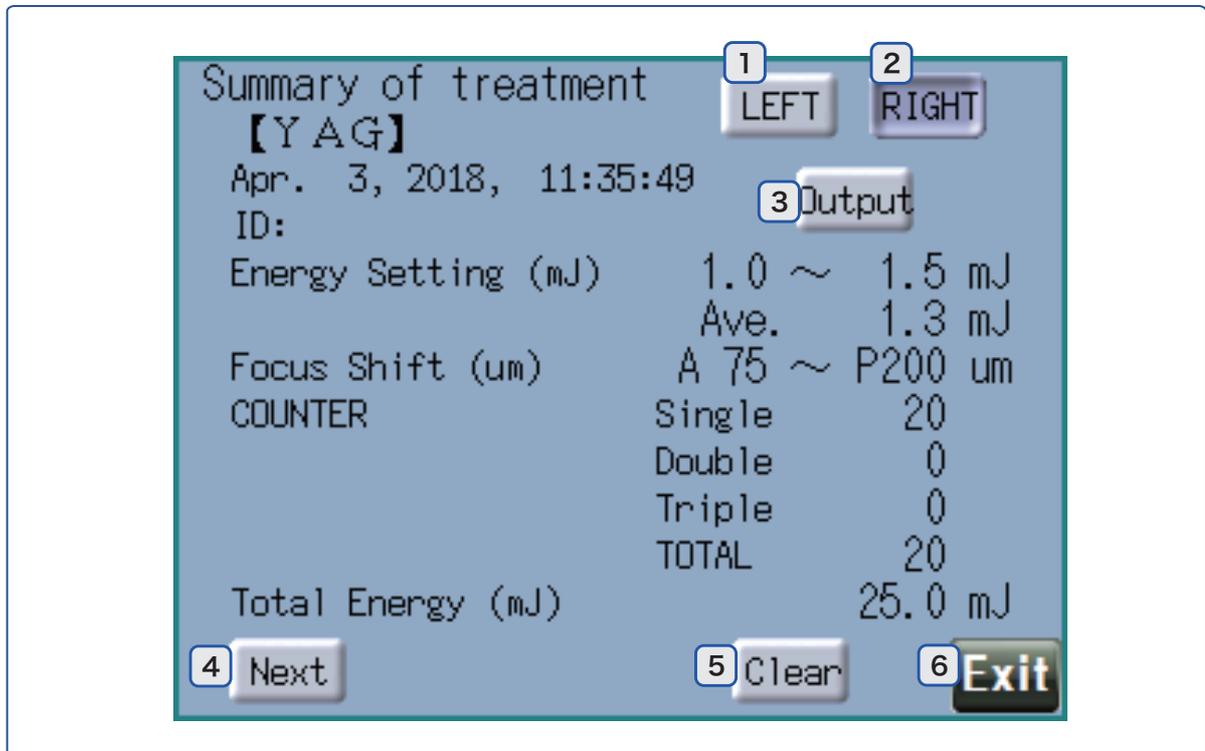

1 LEFT button

Press this button when the treatment beam setting data is for the left eye. When the data is for both eyes, press as well. The button is disabled on the Previous Settings screen.

2 RIGHT button

Press this button when the treatment beam setting data is for the right eye. When the data is for both eyes, press as well. The button is disabled on the Previous Settings screen.

3 Output button

Outputs the treatment beam setting data to the keycard.

It is displayed when or is pressed.

CAUTION

- When both eyes are treated but the summary data needs to be transmitted for each eye, transmit the summary data after the treatment of an eye.
After the treatment of both eyes, the summary data cannot be transmitted for each eye.

4 Next button 

Displays the Previous Settings screen (Summary of treatment).

5 Clear button 

Clears the treatment beam setting and saves it as the Previous Settings.

**Knowledge**


• The following operations save data as the Summary of treatment of the Previous Settings.

- Press .
- Clear the counter on the main screen.
- Shut down the device.

6 Exit button 


Returns to the main screen.

 [“3.4 Changing Device Settings” \(page 50\)](#)

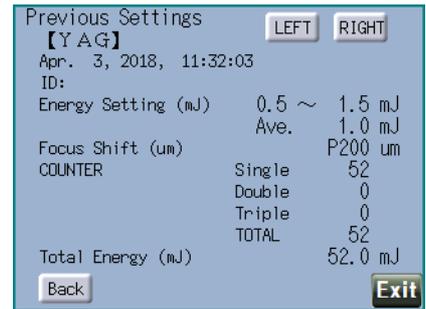**2.11.4 Summary of treatment (SLT mode)**

The Summary of treatment screen is displayed at the time specified in “Summary Display Settings”.

 [“3.4.2 Setting items” \(page 51\)](#)

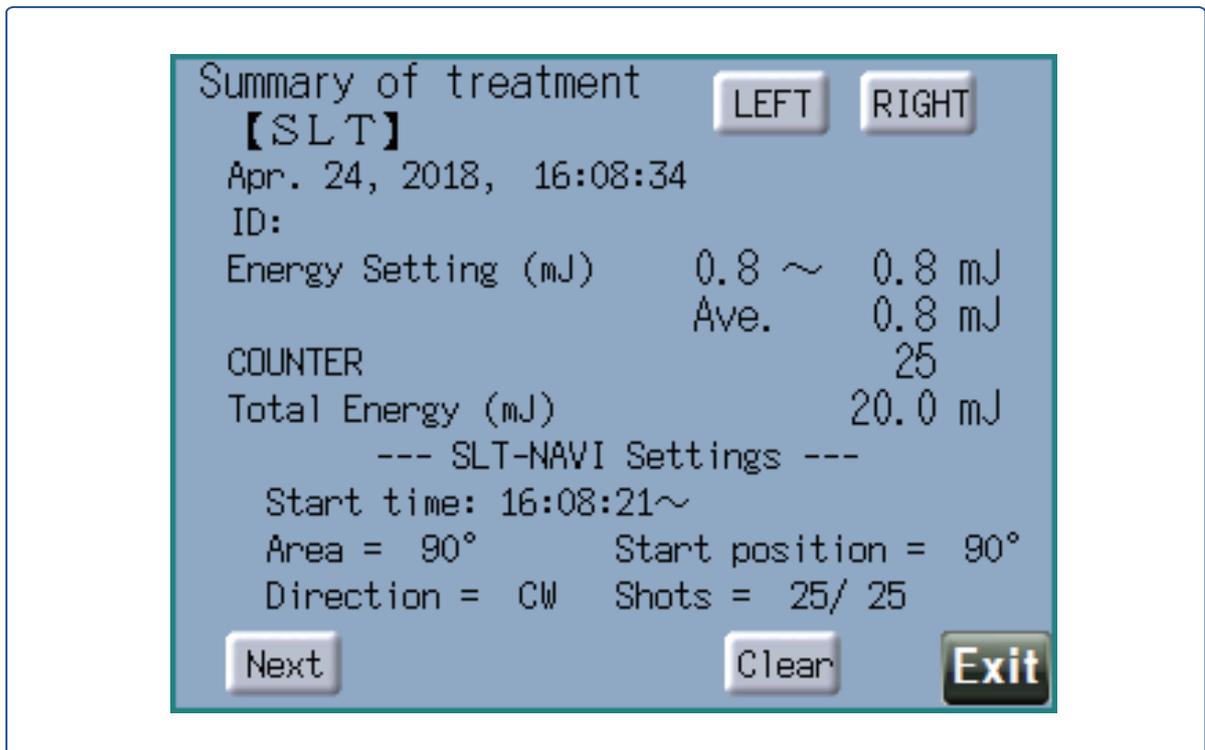
**Knowledge**



- The displayed items on the Summary of treatment screen differ in SLT mode and YAG mode.
- “SLT-NAVI Settings” is not displayed on the screen when SLT-NAVI is inactive.

3

OPERATING PROCEDURE

3.1 Operation Flow

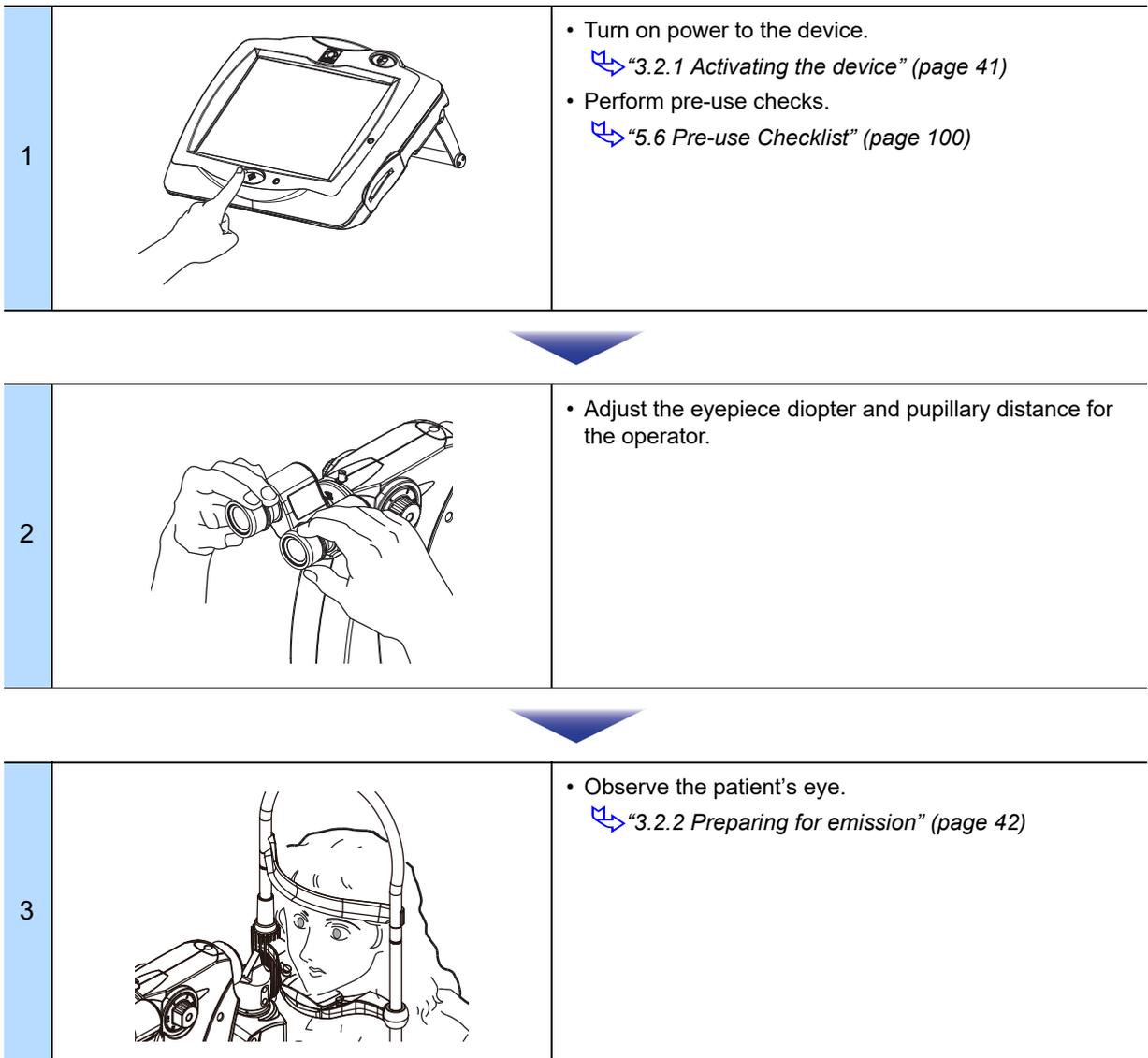

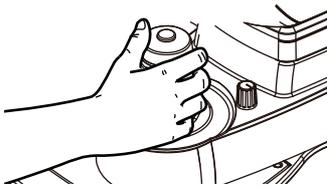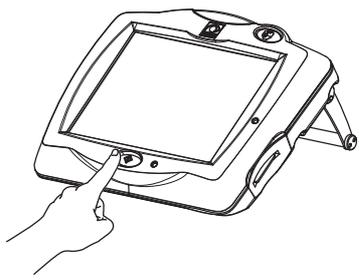4		<ul style="list-style-type: none"> • Set the laser emission conditions.
		
5		<ul style="list-style-type: none"> • Press the status button to set the device to READY mode.
		
6		<ul style="list-style-type: none"> • Emit the treatment beam.
		
7		<ul style="list-style-type: none"> • Turn off power to the device. ↩ "3.2.4 Stopping the device" (page 47)

3.2 Operating Procedure

3.2.1 Activating the device

- 1 Connect the power cord to the power outlet.

WARNING

- Be sure to connect the power plug to a grounded power outlet.
Electric shock or fire may result from device malfunction or electric leakage.


- 2 Have all personnel in the operating room other than the operator and patient wear safety goggles.

WARNING


- Have personnel in the operating room wear safety goggles. Instruct them never to gaze directly at the treatment beam even while wearing the safety goggles.
YAG treatment beam: Wavelength 1,064 nm, OD > 7, 950-1080 D LB6 + IRM LB7 (EN207)
SLT treatment beam: Wavelength 532 nm, OD > 7, 315-532 DIRM LB6 (EN207)

- 3 Activate the device.


- 1) Insert the keycard into the control box.

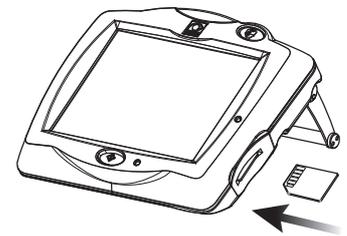

- 2) Turn on (I) the master switch.

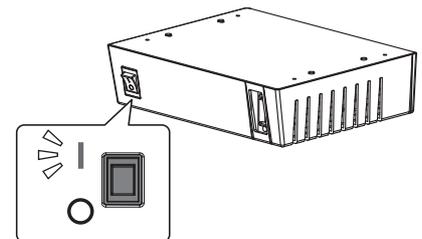

- 3) Press the start button.
The device enters STANDBY mode.

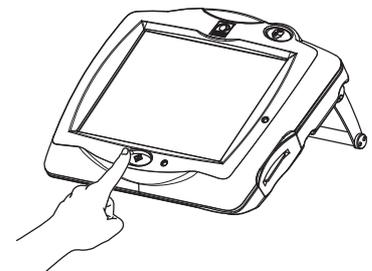

- 4 Perform pre-use checks of the device.


Record the results on the checklist.

 "5.6 Pre-use Checklist" (page 100)

3.2.2 Preparing for emission


1 Adjust the eyepiece diopter and pupillary distance for the operator.

- 1) Remove the cap.

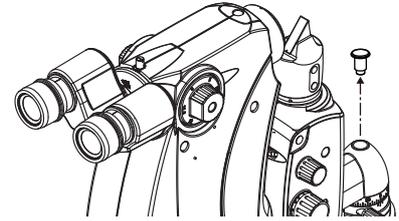


- 2) Insert and turn the focusing rod so that its flat surface faces the microscope.


- 3) Project the illumination light on the focusing rod.

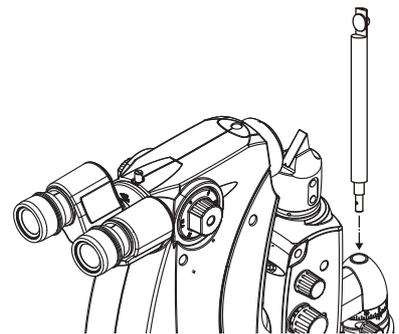

- 4) Fully turn the diopter adjustment rings to the + side and look through the microscope.

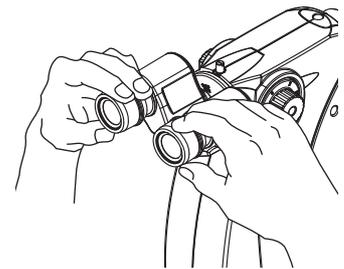

- 5) While observing the slit image with one eye, slowly turn the diopter adjustment ring until the slit image is focused sharply.

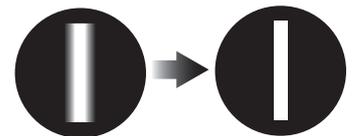


WARNING

- Be sure to adjust the eyepiece diopter for each eye and do not turn the diopter adjustment ring from the - side to the + side.
Inaccurate diopter adjustment may adversely affect laser beam emission.


- 6) In the same manner, adjust the diopter of the other eye.

- 7) Move the binoculars so that the slit image to be observed can be viewed stereoscopically by both eyes.

- 8) Remove the focusing rod and attach the cap to the slit lamp.

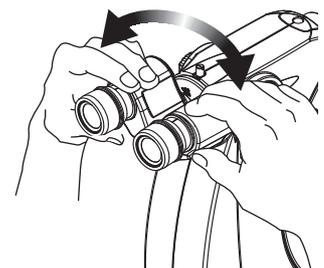

- 2 Press the AIMING button to turn off the aiming beam.

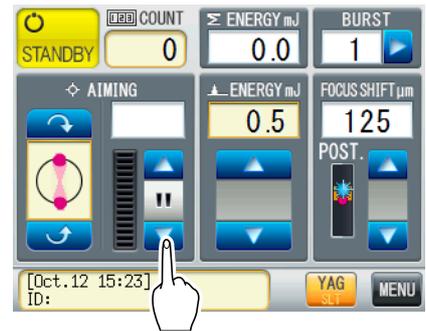


- 3 Conduct patient preparation.


- 1) Pull the joystick all the way toward the operator.

CAUTION


- Be sure to pull the joystick all the way toward the operator before the patient is seated.
This is to prevent the device from coming into contact with the patient's face.

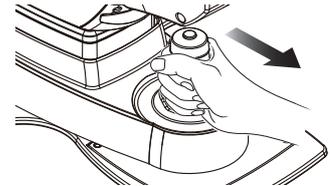


- 2) Instruct the patient to remove their glasses or contact lenses and be seated.
- 3) Instruct the patient to place their chin on the chinrest.
Ask them to rest their forehead on the forehead rest and hold the grips to keep them in a stable posture.

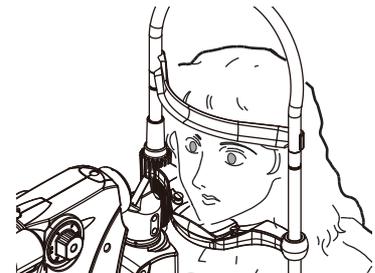

- 4) Rotate the chinrest elevation control to align the level of the patient's eye with the eye level marker.

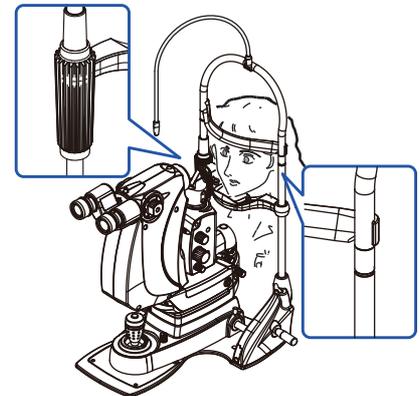

- 5) Fasten the patient's head with the head belt so that the patient's head may not move.
- 6) Instruct the patient to focus on the fixation lamp to stabilize their visual axis.

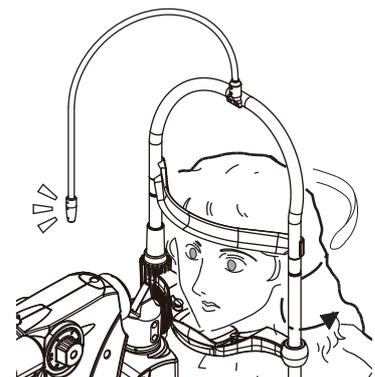


4 Observe the patient's eye.

CAUTION


- Be sure to set the light intensity to the minimum level at the beginning, and raise it as necessary. Be sure to return the light intensity to the minimum level after every examination.
- Minimize the illumination area size as much as possible.
- Maximize the angle between the illumination light and visual axis as much as possible.
- Use an illumination filter as necessary.
- Pay particular attention when projecting the illumination light into the eyes of infants, aphakic patients and patients with eye disease.

1) Manipulate the joystick to project the illumination light on the patient's eye.

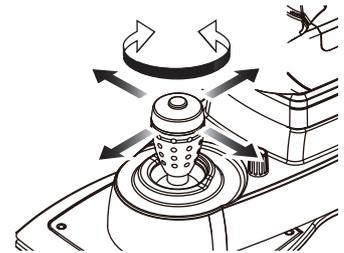

2) Adjust the intensity of the illumination light by turning the illumination control.

3) Manipulate the joystick to focus the slit image on the cornea.

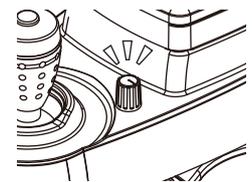


4) Place either a contact lens or hand-held lens in front of the patient's eye.

When using a contact lens, apply topical anesthesia beforehand. If necessary, use a corneal protectant.

5) Look through the microscope to observe the affected areas.


6) As necessary, change the observation conditions.

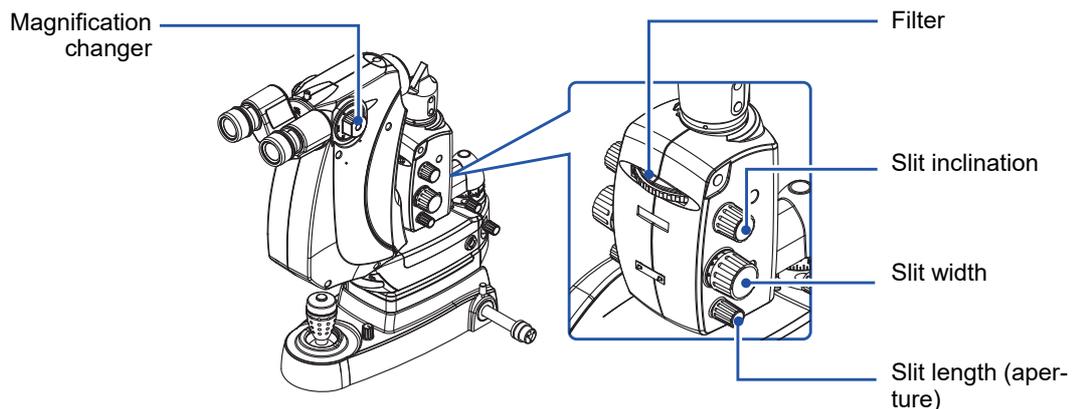

3.2.3 Emitting the treatment beam

1 Set the laser emission conditions.

Changing “ENERGY” or “BURST” automatically performs the test fire.

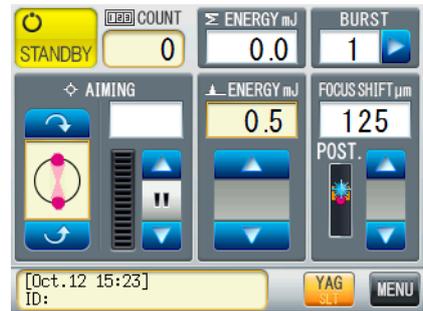

⚠ WARNING

- Check the setting values before emitting the treatment beam. The treatment beam may be emitted in unintended settings.




Knowledge

- During the test fire, a message indicating that no operation is possible appears on the control box.

2 Press the AIMING button to turn on the aiming beam.

⚠ CAUTION


- The aiming beam blinks if the path of the treatment beam is blocked by the illumination tower only when the illumination tower equipped with a tilting function is used. In this case, lower the illumination tower or move the illumination unit arm to the right or left.

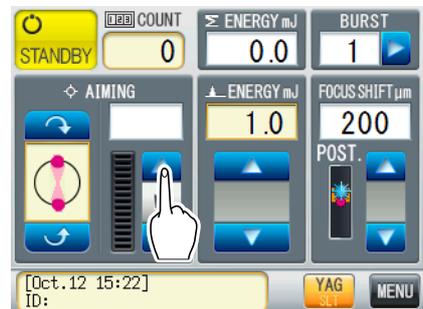

3 Press the status button to set the device to mode.



4 Adjust the joystick and contact lens to align the aiming beam focus with the target position.

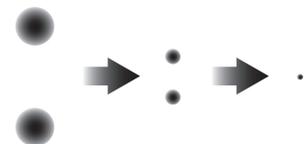

5 Press the hand switch to emit the treatment beam.

⚠ CAUTION

- If the patient's fixation is unstable, emit the treatment beam with great care, paying attention to the position of the aiming beam.




Knowledge

- As a guideline for continuous irradiation, the number of times the hand switch is pressed should not exceed 18 times/minute in Single mode or 6 times/minute in Burst mode. If continuous irradiation that exceeds the above guideline is performed, an error may occur.
- If the laser emission repetition rate exceeds the specified continuous irradiation / repetition time (Single mode: 3 Hz, Burst mode: 1.5 Hz), a message indicating that no operation is possible appears on the control box.

6 When treatment beam emission is complete, press the status button to set the device to mode.




7 Press the AIMING button to turn off the aiming beam.

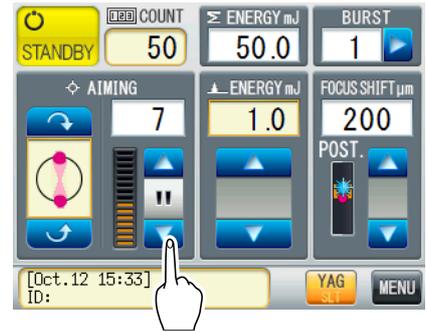

8 Turn off the illumination light by turning the illumination control.


9 Remove the head belt from the patient.

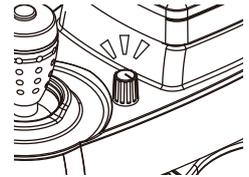

3.2.4 Stopping the device


1 Turn off the device.

- 1) Press the start button.

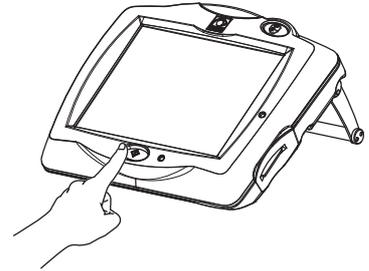


- 2) Press .


- 3) Instruct all personnel to remove their safety goggles.

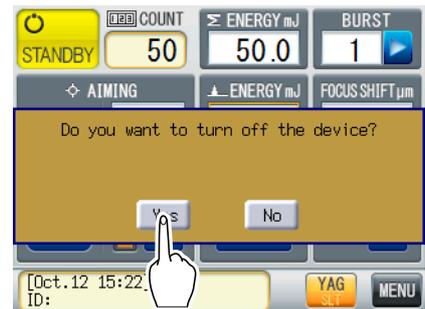

- 4) Turn off () the master switch.

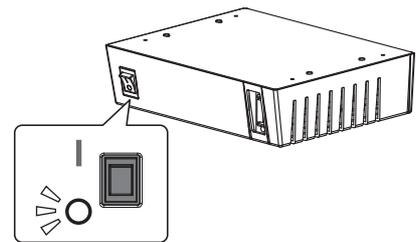

- 5) Remove the keycard and store it in a secure place.

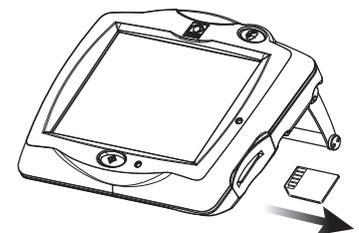

- 2 Disconnect the power cord from the power outlet.

- 3 Clean the forehead rest, chinrest, grips, and head belt.

↳ "4.4.2 Cleaning areas that come into contact with patients" (page 71)

- 4 Place the dust cover on the device.

3.3 Using SLT-NAVI (SLT Mode)


- SLT-NAVI assists the operator in surgery by specifying the laser emission positions and orders before the treatment.




Knowledge

- This function is used to guide the laser emission, not to record the laser emission positions or number of emissions for each area.


1 Set the device to SLT mode.

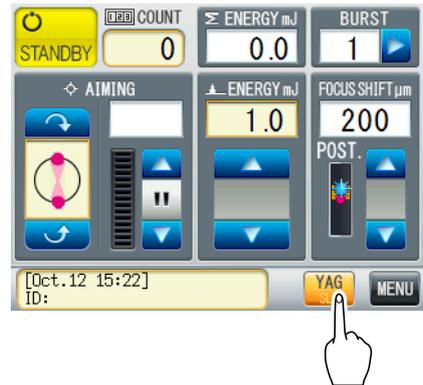


2 Press  .

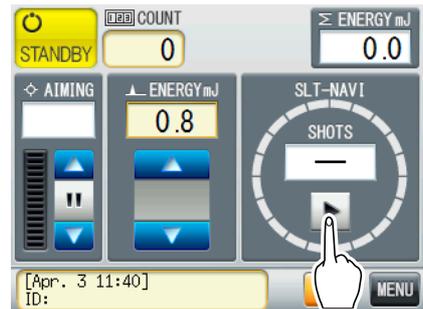


3 Select the emission area.

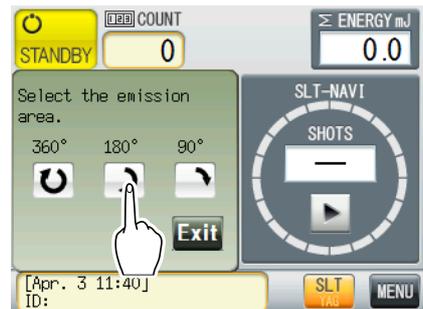

4 Press a light blue area to select the emission start position.

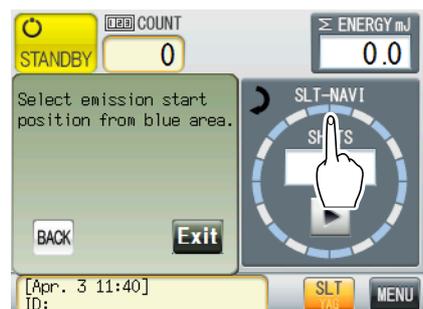

5 Select the emission direction.

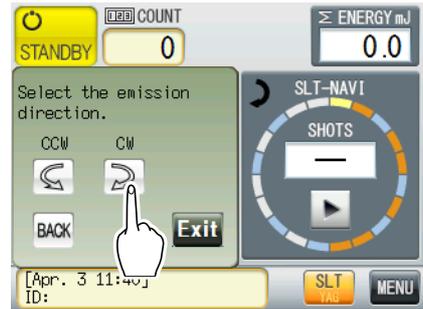

6 Press **OK**.



7 Emit the treatment beam while checking SLT-NAVI.

The procedure for treatment beam emission is the same as that of YAG mode.

↳ "3.2.3 Emitting the treatment beam" (page 45)


Button	Indication/function
Circular shot progress indicator	Orange: Emission target area Blink in yellow and orange: Area that is being irradiated Green: Area that has been irradiated
SHOTS	Number of laser emissions / scheduled number of shots Pressing an area before the treatment beam shot enables changing the scheduled number of shots.
End	End
 	Pause The emission during the pause is not counted.
▶	Restart




3.4 Changing Device Settings

3.4.1 Changing the settings

1 Press **MENU**.

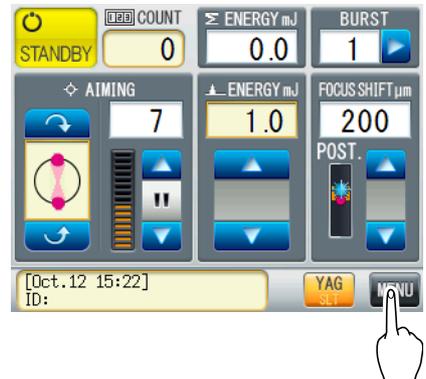


2 Press a button corresponding to the setting item to be changed.

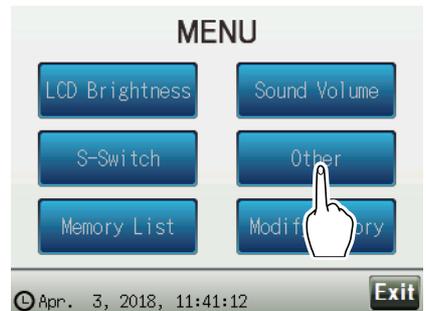

3 Change the setting as desired.

Button	Function
	Selected (cannot be selected if the text is grayed out)
	Returns to the previous screen.
	Goes to the next page.
	Returns to the previous page.

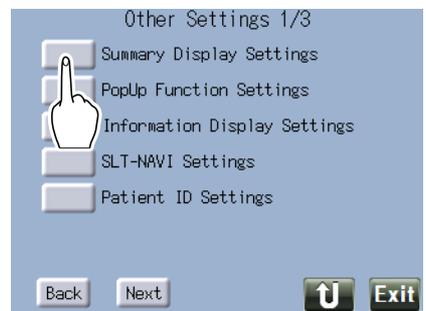

4 Press **Exit**.

The main screen is displayed again. The changed settings are saved.

3.4.2 Setting items



Setting

- Underlined options indicate factory settings.

LCD Brightness		Setting contents
LCD Brightness at STANDBY	<u>1</u> , <u>2</u> , <u>3</u> , <u>4</u> , <u>5</u> , <u>6</u> , <u>7</u>	Select the brightness of the touch screen in STANDBY mode.
LCD Brightness at READY	<u>1</u> , <u>2</u> , <u>3</u> , <u>4</u> , <u>5</u> , <u>6</u> , <u>7</u>	Select the brightness of the touch screen in READY mode.
Sound Volume		Setting contents
Beep Volume:Laser Emission	<u>1</u> , <u>2</u> , <u>3</u>	Select the beep volume when the laser is emitted.
Beep Volume:Button	<u>1</u> , <u>2</u> , <u>3</u>	Select the beep volume when a button is pressed.
S-Switch		Setting contents
Left switch	OFF, READY/STANDBY, AIMING ROTATION START/STOP, BURST , ENERGY UP , ENERGY DOWN , FOCUS SHIFT UP , FOCUS SHIFT DOWN , AIMING POWER UP , AIMING POWER DOWN	Select a function to be assigned to the switch.
Right switch	OFF, READY/STANDBY, AIMING ROTATION START/STOP, BURST , ENERGY UP , ENERGY DOWN , FOCUS SHIFT UP , FOCUS SHIFT DOWN , AIMING POWER UP , AIMING POWER DOWN	Select a function to be assigned to the switch.
Push switch	OFF, <u>READY/STANDBY</u> , AIMING ROTATION START/STOP, BURST , ENERGY UP , ENERGY DOWN , FOCUS SHIFT UP , FOCUS SHIFT DOWN , AIMING POWER UP , AIMING POWER DOWN	Select a function to be assigned to the switch.
Default switch		Restore the functions assigned to the switch to the default settings.
Other		Setting contents
Summary Display Settings	<u>Not Display</u> , STANDBY , COUNTER Reset	Select when to display the Summary of treatment screen.
PopUP Function Settings	No , <u>Yes</u>	Select whether to display the pop-up window for selection of the aiming beam angle and energy value.

Other

Setting contents

Information Display Settings

Line 1: **Not Display, Patient ID, Date**
 Line 2: **Not Display, Patient ID, Date, Memory Name**

Select the items to be displayed in the information display field.

SLT-NAVI Settings (SLT mode)

90, 100, 110, 120, 130, 140, 150, 160

Select the number of shots for entire circumference when SLT-NAVI is active.

Patient ID Settings

Enter a patient ID with the on-screen keyboard.

Confirm it by pressing **Enter** .

Pressing a button other than “a to z”, “A to Z”, “0 to 9”, and “-” displays “~” (tilde).



Data Input settings

Cannot be selected.

Data Output settings

Cannot be selected.


AIMING Rotation Increment Settings (YAG mode)

Step1, Step2, Step3

Select the increments for rotating the aiming beam.



displays the increments for rotation.

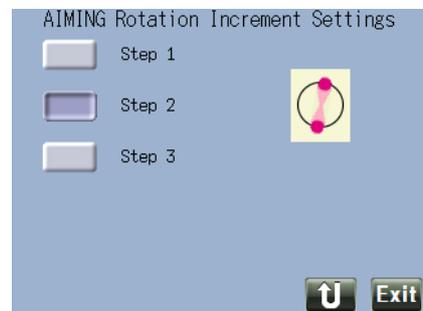

Σ ENERGY Display Settings

No, Yes


Select whether to display the integrated energy value.

Clock Settings

Sets the date and time.

Enter values in the numeric keypad displayed when a button is pressed.

Values from 2017 to 2099 can be entered as a year.

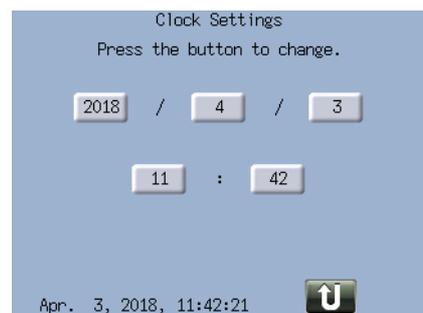

Other

Setting contents

Display Language Settings

English, Japanese

Select the display language.

Trigger Switch Selection


Hand Switch, Foot Switch (optional)

Select the trigger switch used to emit the laser.

Setting Range Limit (Selected at Start-up) (YAG mode)

No, Yes


Select whether to enable Limit mode at device start-up.

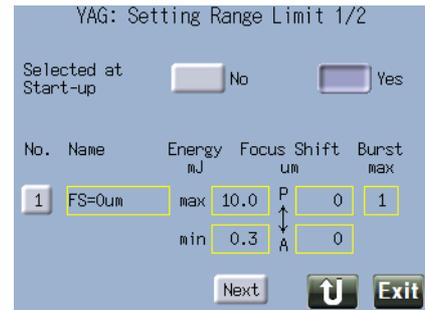

Setting Range Limit (Name) (YAG mode)

Press a yellow frame to enter a name.

Confirm it by pressing **Enter**.

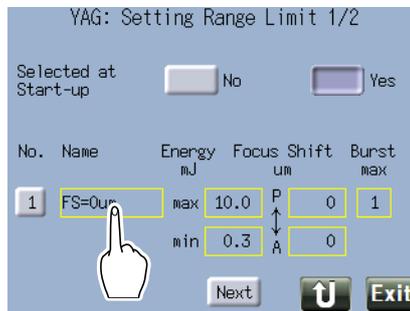

Other

Setting contents

**Setting Range Limit
(Energy, Focus Shift,
Burst)
(YAG mode)**

Press a yellow frame and enter a value with the numeric keypad.



Energy: Set in the range from 0.3 (mJ) to 10.0 (mJ).

Focus Shift: Set in the range from P500 (Post: 500 um) to A500 (Ant: 500 um). Values can be entered in 25 um increments.

Burst: Set in the range from 1 to 3. It can be set in combination with the following Energy.

- Burst 1: 0.3 to 10.0 mJ
- Burst 2: 0.3 to 9.0 mJ
- Burst 3: 0.3 to 8.0 mJ

Confirm it by pressing **Enter**.

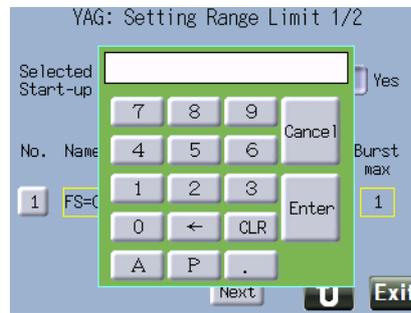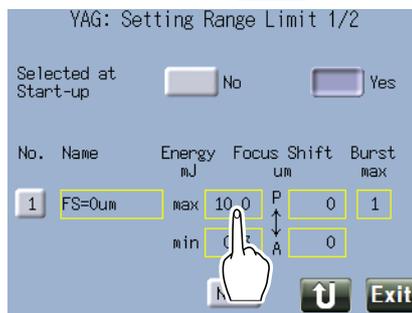

CALIBRATION

Calibrate the laser energy output.

3.5 Using Limit Mode (YAG Mode)

- In Limit mode, the setting ranges of Energy, Focus Shift, and Burst are limited to use the device more safely. Up to three settings can be registered to Limit mode.


3.5.1 Using Limit mode from start-up



Setting


- Set "Selected at Start-up" of "Setting Range Limit" to **Yes** in advance.

↳ "3.4 Changing Device Settings" (page 50)

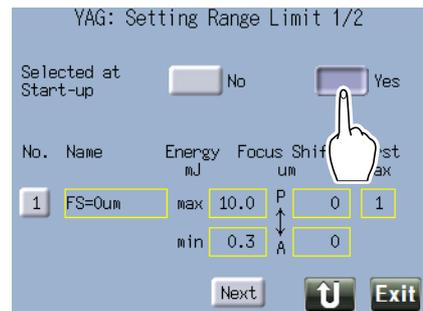

1 Activate the device.

↳ "3.2.1 Activating the device" (page 41)

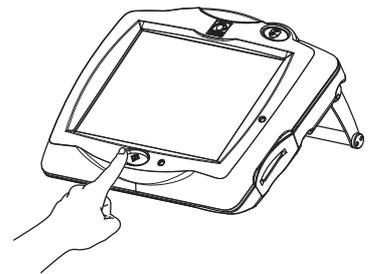

2 Press **1**, **2**, or **3**.

Button	Function
1 to 3	Activates Limit mode.
SLT	Activates normal SLT mode.
Exit	Activates normal YAG mode.




3 The Limit mode is activated.

"Limit Mode" is displayed in the information display field.




3.5.2 Using Limit mode after activating YAG mode


1 Press **MENU**.

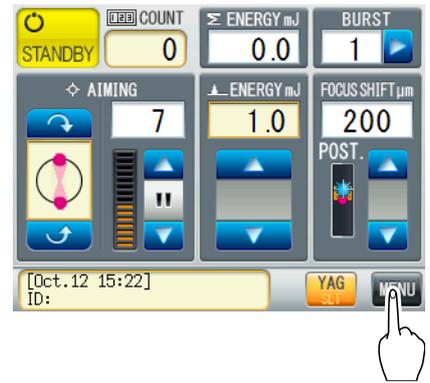

2 Press "Other".

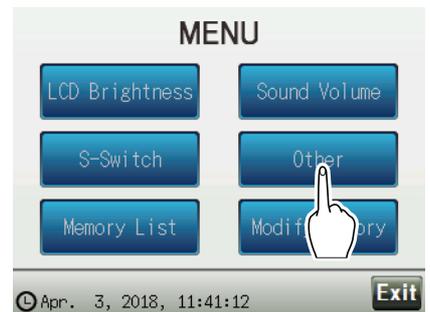


3 Press **Next** to display the Other Settings 3/3 screen.

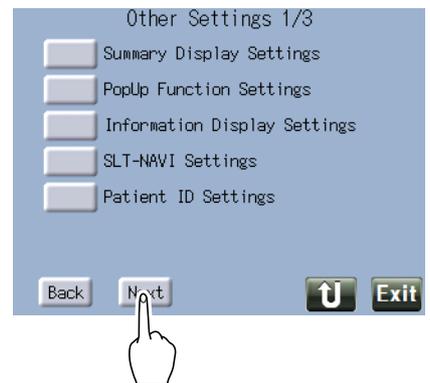

4 Press "Setting Range Limit".




5 Press the No. button of the desired Limit mode.

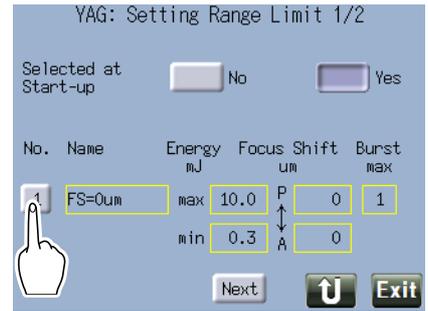

6 Press **ON**.

Button	Function
ON	Turns on Limit mode.
OFF	Turns off Limit mode.
Clear	Deletes the data. Energy, Focus Shift, and Burst are changed to the default value.
Exit	Closes the pop-up window.

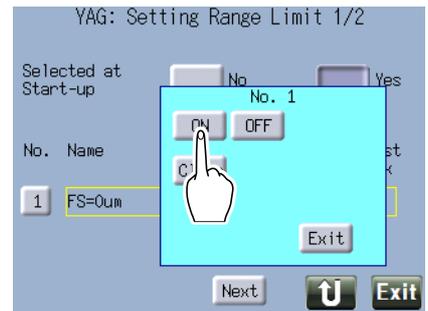

7 Press **Exit**.


8 The Limit mode is activated.

“Limit Mode” is displayed in the information display field.




3.5.3 Deactivating Limit mode

1 Press the information display field.

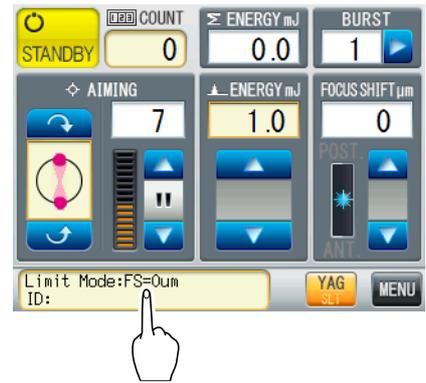

2 Press .

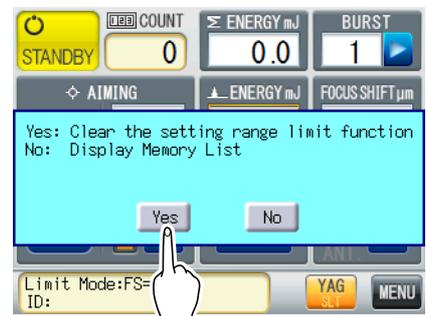


3.6 Using Memory


- Treatment beam settings frequently used can be registered to the memory for easy setting.

3.6.1 Loading memory


1 Press **MENU**.

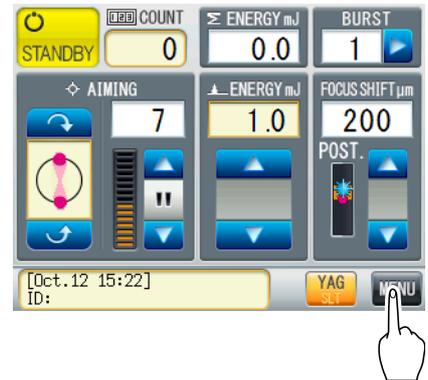

2 Press "Memory List".

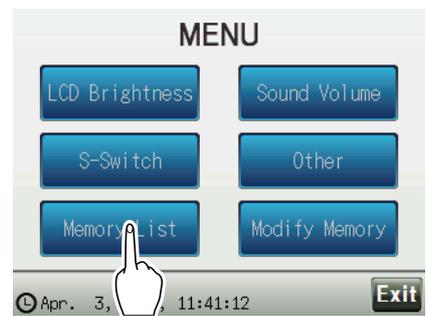


3 Press the No. button of the memory to be loaded.

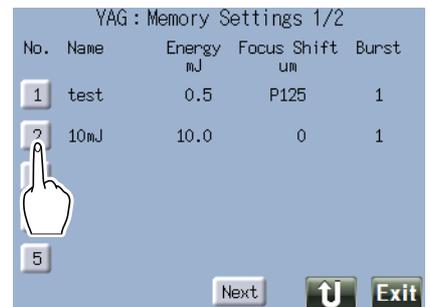

Tip


- Pressing the information display field displays the memory list in a single action.
- When Limit mode is activated, select the memory list or deactivate Limit mode.

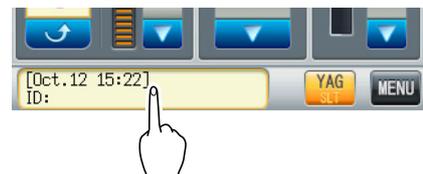


3.6.2 Registering memory

The treatment beam setting that is displayed on the main screen can be registered to memory. Registered memory can be loaded from the memory list.


1 Press **MENU**.

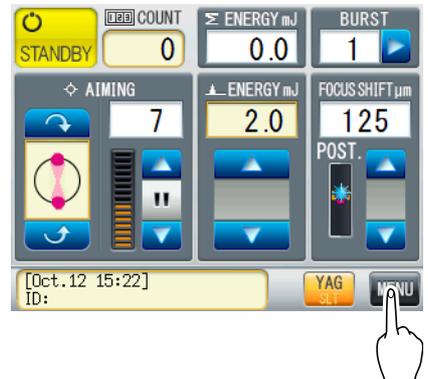


2 Press "Modify Memory".

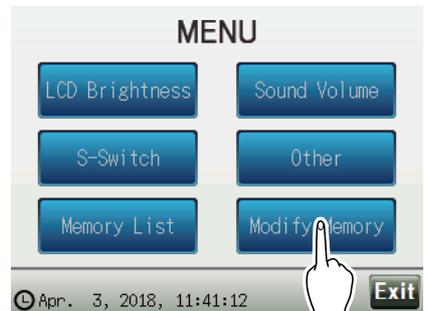


3 Press the No. button of the memory to be registered.

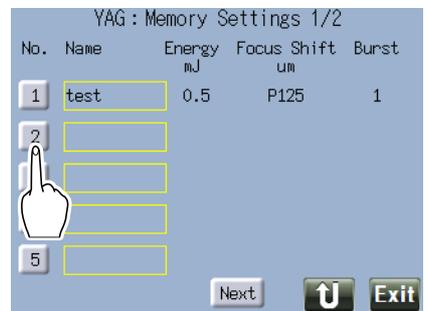


4 Press **Write**.

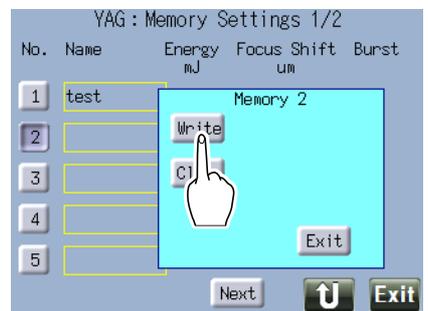

5 Press **Yes** .

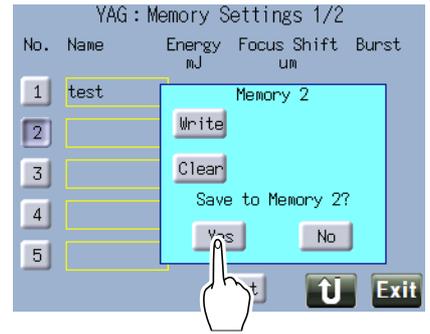

6 Press the yellow frame of the memory to be registered.

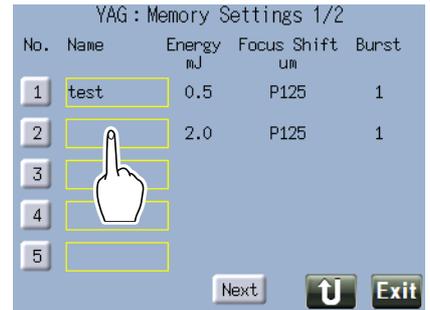


3


7 Enter the name.



8 Press **Enter** .

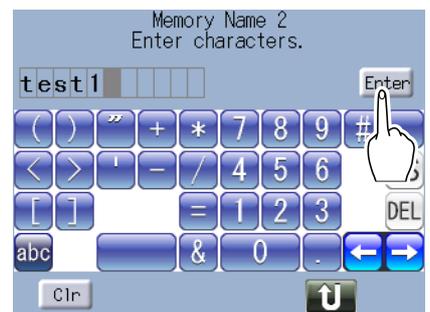


9 Confirm that the treatment beam setting data is registered.

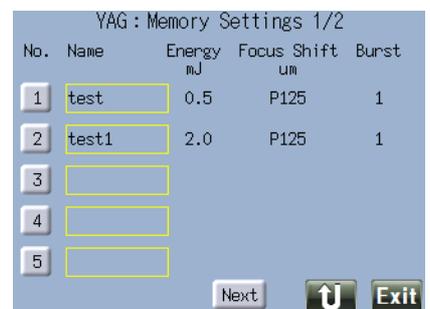


3.6.3 Deleting memory


1 Press **MENU** .

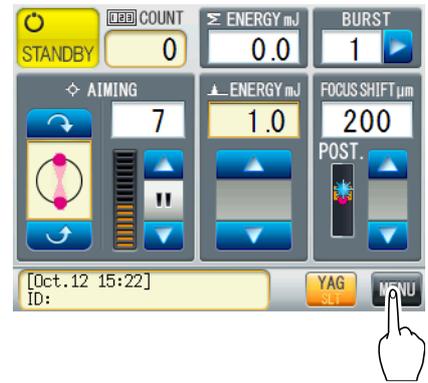


2 Press “Modify Memory”.

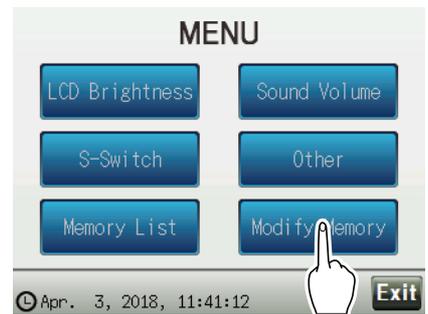


3 Press the No. button of the memory to be deleted.

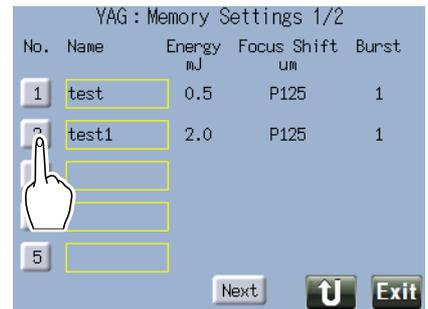

4 Press **Clear** .

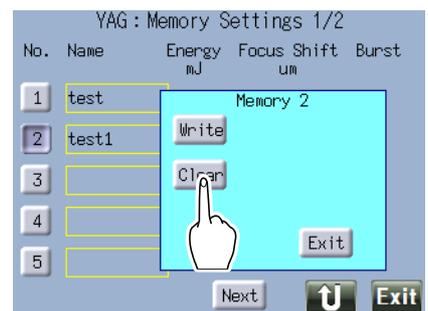

5 Press **Yes** .

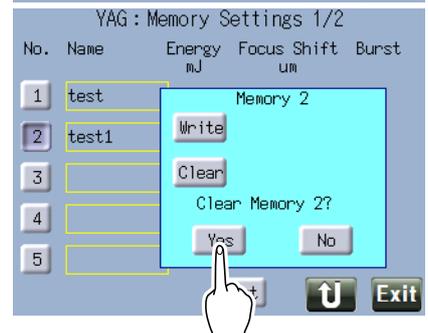

- 6** Confirm that the treatment beam setting data is deleted.

YAG: Memory Settings 1/2

No.	Name	Energy mJ	Focus Shift um	Burst
1	test	0.5	P125	1
2				
3				
4				
5				

Next  Exit

4


MAINTENANCE

4.1 Troubleshooting

If the device does not function properly, attempt to correct the problem according to the following table before contacting NIDEK or your authorized distributor.

If the symptom cannot be corrected by the actions shown in the following table, contact NIDEK or your authorized distributor.

When	Remedy
The power indicator does not illuminate even though the master switch is on (I).	<ul style="list-style-type: none"> The power cord may not be correctly connected. Check the power cord. Check whether proper voltage is applied to the power outlet.
The screen of the control box blacked out.	<ul style="list-style-type: none"> The device may be in sleep mode. In sleep mode, the status indicator slowly blinks in orange. Touch the touch screen or press the start button of the control box.
A message "The KEY card is not inserted. Please enter password." is displayed.	<ul style="list-style-type: none"> The keycard is not inserted. Insert the keycard.
The illumination light is not projected.	<ul style="list-style-type: none"> The illumination light may be turned off. Adjust the intensity of the illumination light by turning the illumination control. The slit may be closed. Increase the slit width by turning the slit width control. It may be blocked by the filter. Adjust the filter position by turning the filter changer. The cable may not be correctly connected to the connector box. Check the cable.
The fixation lamp does not illuminate.	<ul style="list-style-type: none"> The fixation lamp cord may not be correctly connected to the connector box. Check the fixation lamp cord.
The device does not move up or down.	<ul style="list-style-type: none"> The device has reached the upper or lower limit of its movement range. Rotate the joystick in an opposite direction.
The aiming beam is not projected.	<ul style="list-style-type: none"> The aiming beam may be turned off. Press the AIMING button to turn on the aiming beam. The aiming beam may be blocked by the illumination tower. Move the illumination unit arm to the right or left to adjust the position so that the aiming beam is not blocked. In another way, lower the illumination tower to the position where the aiming beam is not blocked (when the illumination tower equipped with a tilting function is used).
The aiming beam blinks.	<ul style="list-style-type: none"> The aiming beam blinks if the treatment beam is blocked due to the position of the illumination tower (when the illumination tower equipped with a tilting function is used). Move the illumination unit arm to the right or left to adjust the position so that the treatment beam is not blocked. In another way, lower the illumination tower to the lowest level.

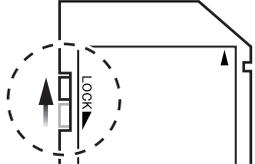When	Remedy
<p>When the illumination tower is replaced with the optional one or that of the combination delivery unit, the laser beam is blocked by the illumination tower.</p>	<ul style="list-style-type: none"> • Move the illumination unit arm to the right and left to adjust the position so that the laser beam is not blocked.
<p>The treatment beam cannot be emitted.</p>	<ul style="list-style-type: none"> • The treatment beam is not emitted if it is blocked due to the position of the illumination tower. Move the illumination unit arm to the right or left to adjust the position so that the treatment beam is not blocked. In another way, lower the illumination tower to the lowest level (when the illumination tower equipped with a tilting function is used).
<p>Writing to the keycard is not possible.</p>	<ul style="list-style-type: none"> • The write-protection switch may be in the LOCK position. Release write-protection. 

4.2 Error Messages and Remedies

If one of the following error codes is displayed on the screen, follow the instructions of the message.

Notify NIDEK of the error code, message number, and serial number of your device so that NIDEK can offer appropriate service.

◆ Warning message

Message No.	Message
• Related to control box	
601	CB EEPROM: Check Sum error CB EEPROM was initialized.
631	Voltage reduction for clock data backup Operate the device for more than 5 minutes for battery charge.
633	Clock data initialized Set the clock.
641	KEY card: Error KEY card has fault. Recording function is suspended. Do not use this KEY card.
644	KEY card: Insufficient free space Free up space of this KEY card.


◆ Error Messages

The following are errors related to the mechanism inside the device. Turn off the device, then restart the device. If the problem persists, contact NIDEK or your authorized distributor.			
• Errors related to main body			
ERR1	High temperature in the system	ERR5	Malfunction of YAG shutter
ERR6	Malfunction of SLT shutter	ERR7	Malfunction of safety shutter
ERR10	Insufficient laser emission	ERR12	Excessive laser emission
ERR13	Incorrect laser emission	ERR14	Laser emission failed
ERR15	Unstable YAG laser output	ERR16	Unstable SLT laser output
ERR30	Laser power supply malfunction	ERR31	Improper PFN voltage of laser power supply
ERR32	Improper PFN voltage for automatic discharge	ERR50	Energy monitor malfunction
ERR51	Unable to set YAG energy	ERR52	Unable to set focus shift
ERR53	Trigger switch input signal failure	ERR54	Pulse number setting signal failure
ERR56	Unable to set aiming beam rotation	ERR61	Unable to switch between YAG and SLT
ERR62	Unable to set SLT energy	ERR65	Unable to control KTP temperature
ERR71	Improper reference voltage	ERR80	Slit lamp malfunction
ERR90	Program data failure	ERR91	EEPROM data failure
ERR93	RAM data failure	ERR94	KEYLOCK error
ERR95	CPU communication error	ERR131	PFN voltage upper limit error of laser power

ERR165	KTP temp control upper/lower limit error	ERR191	EEPROM data editing failure
ERR192	System reset error		
• Errors related to control box			
ERR602	EEPROM:Access error	ERR612	Flash memory-Appli:Check Sum error
ERR613	Process:Program Check Sum error	ERR614	Process:Data mismatch error
ERR615	Process:Calculation error	ERR616	Process:Sequence error
ERR622	KEY card :Use of unauthorized KEY card.	ERR632	RTC:Access error
ERR642	KEY card :KEY card is not inserted.	ERR643	KEY card :Write-protected.
ERR651	No response error.	ERR652	Communication error.
ERR661	Cannot turn off power.	ERR671	LCD BL driver:Error
ERR681	VerUP:File Error	ERR683	VerUP CB:Erase Error
ERR684	VerUP CB:Write Error	ERR685	VerUP CB:Verify Error
ERR692	VerUP MAIN:Bootmode Error	ERR693	VerUP MAIN:Erase Error
ERR694	VerUP MAIN:Write Error	ERR695	VerUP MAIN:Verify Error

4.3 Replacing Consumables

4.3.1 Consumable list

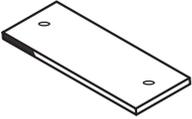
Part name (part number)	Appearance	Remarks
Chinrest paper (32903-M047)		1 pack
Fuse For 100 V regions: 80402-02148 For 200 V regions: 80402-02043		For optional optical table For 100 V regions: 125 V T 10 A For 200 V regions: 250 V T 5 A
Slide plate (YC020-M216)		For replacement, contact NIDEK or your authorized distributor.

- After replacing any consumables, be sure to restock them with spares.

4.3.2 Attaching chinrest paper



Attach an appropriate amount of chinrest paper onto the chinrest using the fixing pins.

Take out a stack of chinrest paper from the chinrest paper pack so that the thickness of the stack is 6 mm or less, then attach the stack of the chinrest paper.

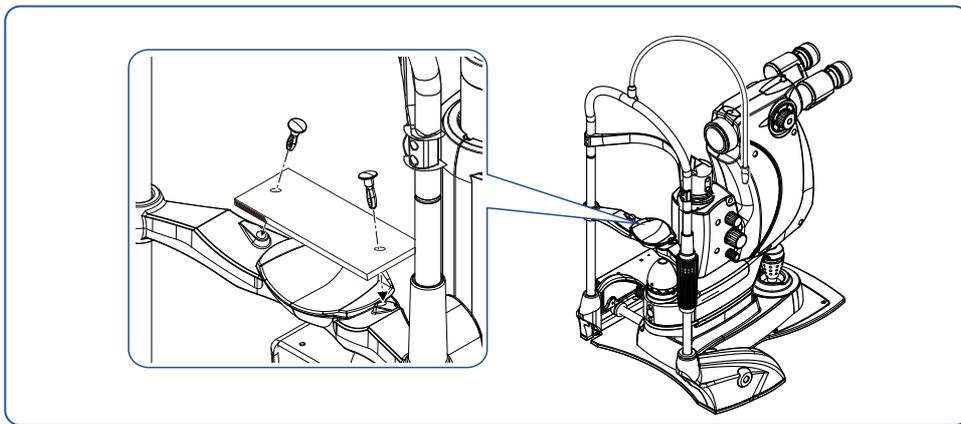

4.3.3 Replacing fuses of optional optical table

- 1** Turn off () the master switch.
- 2** Disconnect the power cord from the power outlet.
- 3** Remove the fuse holders.
Turn the fuse holders counterclockwise with a flatblade screwdriver while applying pressure.
- 4** Replace the fuses.
Replace the two fuses at the same time.
- 5** Reassemble the parts in the reverse order.

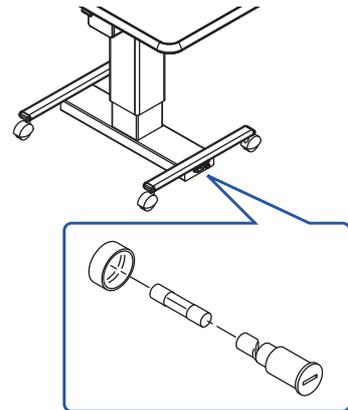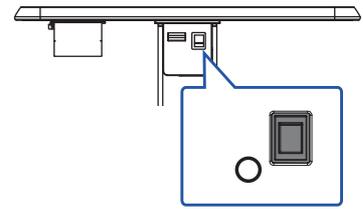

4.4 Cleaning

4.4.1 Cleaning the device exterior

CAUTION

- Disconnect the power cord from the power outlet.
This is to prevent electric shock.
- Never use an organic solvent (such as paint thinner and benzine) or abrasive cleanser.
It may damage the surface of the device.
- Gently wipe the touch screen surface with a soft cloth.
It may scratch the touch screen. Device malfunction may also result.
- Never use an overly wet sponge or cloth.
Water may leak into the interior of the device resulting in malfunction.

1 Turn the power off ().

2 Disconnect the power cord from the power outlet.

3 Wipe off the stains with a soft cloth.


For persistent stains, soak the cloth in a neutral detergent, wring well, and wipe. Then wipe them off with a soft, dry cloth.

4 Let it dry in a well-ventilated area.

4.4.2 Cleaning areas that come into contact with patients

Clean the forehead rest, chinrest, grips, and head belt for each patient. Use alcoholic disinfectants if necessary.

When using the chinrest paper, remove one sheet of paper.

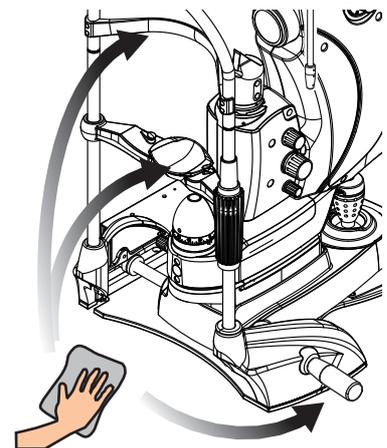

4.4.3 Cleaning optical parts

 **CAUTION**

- Wipe up to clean the optical parts.
It the optical parts are not clean, the laser energy output may be unstable.
 - Take care not to scratch or smear the optical parts, or accumulate dust on them.
This prevents reduction in laser performance.
-
-

- 1** Remove any dust by using the blower.
- 2** Wrap lens cleaning paper around a thin stick (or use cotton swab / gauze), damp it with alcohol, and wipe the optical parts.

Gently wipe the lens circularly from the center to periphery.



Check

- Do not use a metal stick or something made of hard materials.

- 3** Confirm that there is no dust or smudge.

If the optical parts are not cleaned satisfactorily, repeat the cleaning procedure with new lens cleaning paper until they become clean.

4.5 Checking Optical Axis of Treatment Beam (YAG mode)


- Confirm that the aiming beam and treatment beam are coaxial.

1 Prepare the following things.

- Test paper (business card-sized black paper)
- Cellophane tape


2 Activate the device.

↳ “3.2.1 Activating the device” (page 41)

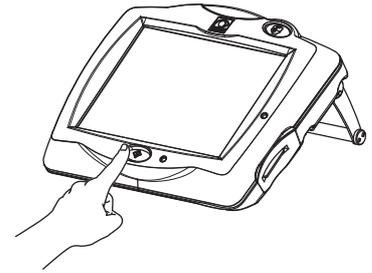

3 Adjust the eyepiece diopter and pupillary distance for the operator.

↳ “3.2.2 Preparing for emission” (page 42)

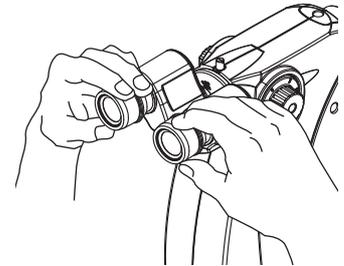


4 Specify the settings of the device.


Item	Setting
Microscope magnification	32x
ENERGY	0.3 to 0.5 mJ
FOCUS SHIFT	0 μm

5 Attach the test paper to the forehead rest with cellophane tape.


Position it perpendicularly to the optical axis as much as possible.

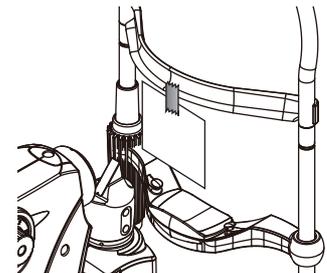

6 Emit the aiming beam onto the test paper to adjust the focus.

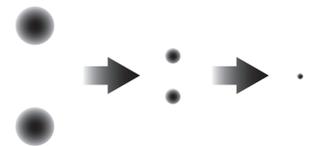

7 Press the status button to set the device to mode.

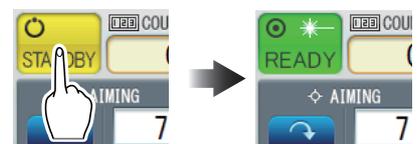

- 8** Emit the treatment beam.
- 9** Confirm that the aiming beam is aligned with the burn on the test paper.



4.6 Checking Optical Axis of Treatment Beam (SLT mode)


- Confirm that the aiming beam and treatment beam are coaxial.

1 Prepare the following things.

- Test paper (business card-sized black paper)
- Cellophane tape


2 Activate the device.

→ “3.2.1 Activating the device” (page 41)

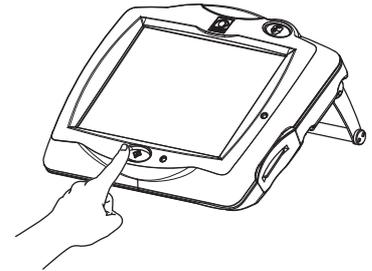

3 Adjust the eyepiece diopter and pupillary distance for the operator.

→ “3.2.2 Preparing for emission” (page 42)

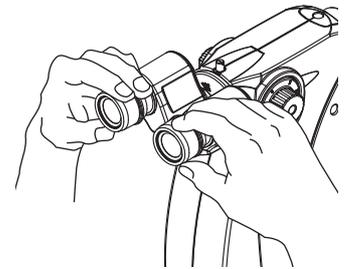


4 Specify the settings of the device.

Item	Setting
Microscope magnification	32x
ENERGY	0.3 to 0.5 mJ

5 Attach the test paper to the forehead rest with cellophane tape.

Position it perpendicularly to the optical axis as much as possible.

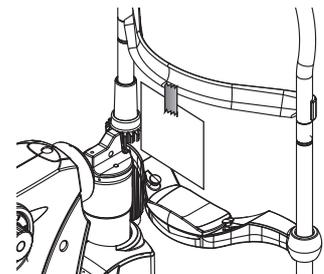

6 Emit the aiming beam onto the test paper to adjust the focus.

7 Check the aiming beam.

- The intensity is even.
- The intensity is not lowered. The periphery is not obscured.
- The outline of the aiming beam is clear at the focal position.



8 Press the status button to set the device to mode.



- 9** Emit the treatment beam.
- 10** Confirm that the aiming beam is aligned with the burn on the test paper.



4.7 Measuring and Calibrating Treatment Beam Energy Output

- If the energy output of the treatment beam is out of tolerance, calibration is required.

An energy meter is necessary for energy output measurement and calibration.


CAUTION

- Only service personnel trained by NIDEK or NIDEK distributor are allowed to calibrate the laser energy output.
- Follow the procedure at least once a year.

4.7.1 Measuring treatment beam energy output


1 Mount an energy meter.

- 1) Secure the detector of the energy meter to the chinrest.
Align the height of the center of the detector's receiving surface with the eye level marker.
- 2) Connect the detector cable to the energy meter.

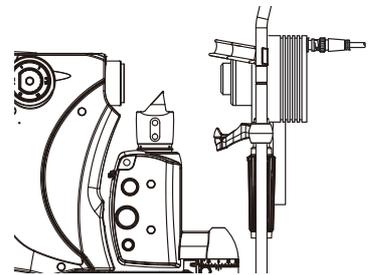

2 Activate the device.

 "3.2.1 Activating the device" (page 41)

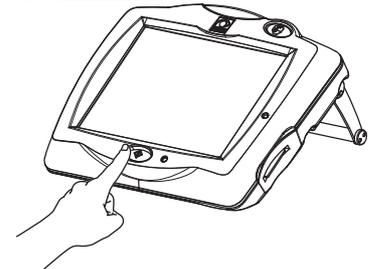


3 Align the aiming beam with the receiving surface of the detector.

- 1) Turn off the illumination light by turning the illumination control.
- 2) Project the aiming beam on the receiving surface of the detector.



3) Manipulate the joystick to move the objective lens of the device from 80 to 90 mm^{*1} away from the receiving surface of the detector.

If the laser beam is focused on the receiving surface, the treatment beam may damage the receiving surface.

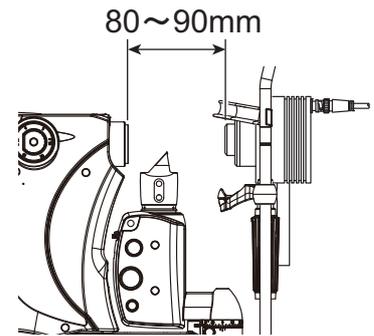

4 Measure the treatment beam for all the combinations in the table.

Mode	ENERGY	BURST	Tolerance
YAG mode	0.3 mJ	1, 2, 3	±20%. For energy output stability, see "5.1 Specifications" (page 87).
	0.5 mJ	1, 2, 3	
	1.0 mJ	1, 2, 3	
	3.0 mJ	1, 2, 3	
	5.0 mJ	1, 2, 3	
	7.0 mJ	1, 2, 3	
	8.0 mJ	3	
	9.0 mJ	2	
SLT mode	10.0 mJ	1	±20%. For energy output stability, see "5.1 Specifications" (page 87).
	0.3 mJ	-	
	0.5 mJ	-	
	1.0 mJ	-	
	2.0 mJ	-	
	3.0 mJ	-	

5 Record the measurement values.

↳ "5.5 Laser Beam Output Calibration Record" (page 98)

6 When the measurement values are out of the tolerance, the energy output needs to be calibrated.


↳ "4.7.2 Calibrating treatment beam energy output" (page 79)

*1. Leave this space when using an energy meter PE50-DIF-C manufactured by OPHIR. The required distance depends on energy meters. Decide the distance necessary for the energy meter to be used.

4.7.2 Calibrating treatment beam energy output

- 1 Change the mode for the treatment beam to be calibrated. (YC-200 S plus)

The treatment beam to be calibrated corresponds to the mode when CALIBRATION mode is activated.

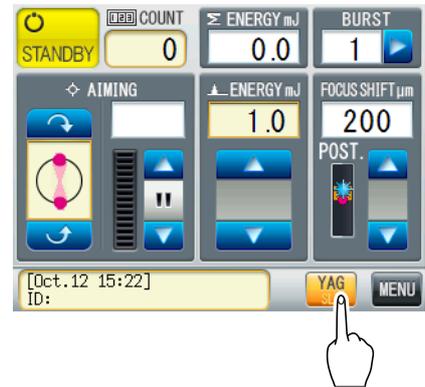

- 2 Press the AIMING button to turn on the aiming beam.




- 3 Activate CALIBRATION mode.

[↩ "3.4 Changing Device Settings" \(page 50\)](#)

Only administrators qualified by NIDEK are informed of the password.




- 4 Press "Energy".

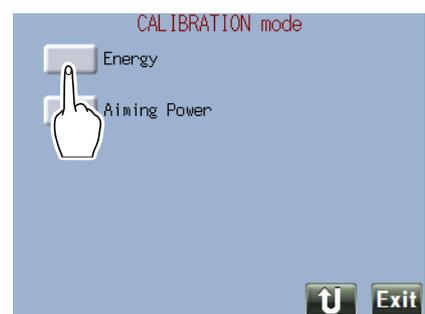


- 5 Change the parameters for 0.3 mJ energy value.

1) Press "0.3 mJ".

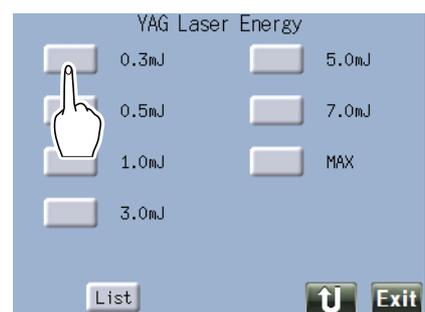

2) Press “STANDBY”.

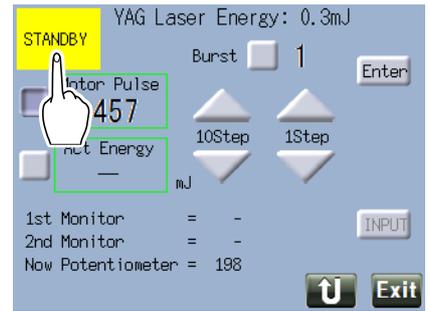

3) Emit the treatment beam and measure the energy.




4) Confirm that the measurement value is within 0.3 mJ \pm 0.01 mJ.

- The measurement value is within 0.3 mJ \pm 0.01 mJ. Proceed to *Step 5*).
- The measurement value is not within 0.3 mJ \pm 0.01 mJ.

Change the parameter for “Motor Pulse” and press **Enter**.

Repeat from *Step 3*) to *Step 4*) until the value is within the criteria.

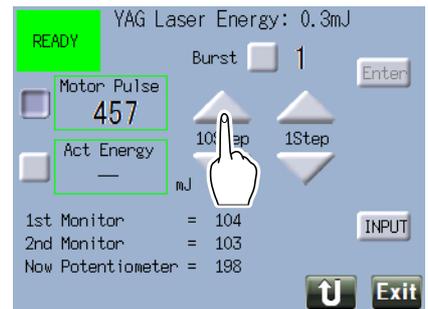


5) Press [Input] to save the parameter for “Motor Pulse”



Knowledge

- [Input] can be pressed after the laser is emitted.

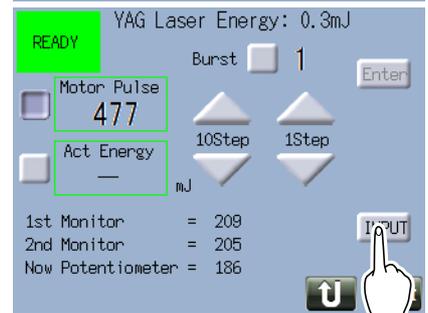

6) Change “Burst” to 2.

Repeat the treatment beam emission and change the parameter for “Motor Pulse” until the measurement value is within 0.6 mJ \pm 0.02 mJ.

When the measurement value is within the criteria, press [Input] to save the parameter.



7) Change "Burst" to 3.

Repeat the treatment beam emission and change the parameter for "Motor Pulse" until the measurement value is within 0.9 mJ \pm 0.03 mJ.

When the measurement value is within the criteria, press [Input] to save the parameter.




8) Press .




Knowledge

- When the parameter for 0.3 mJ is changed, "Act Energy" cannot be selected.

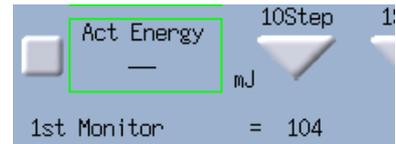


6 Change the parameters for 0.5 mJ energy value.

1) Press "0.5 mJ".

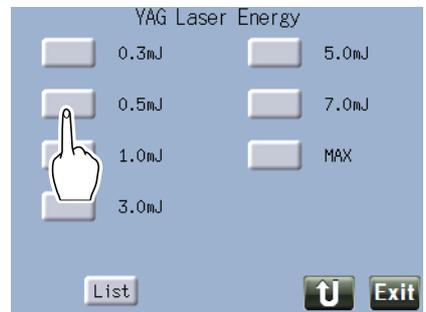

2) Press "STANDBY".

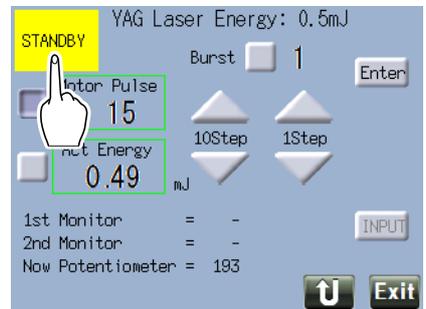

3) Emit the treatment beam and measure the energy.




4) Confirm that the measurement value is within $\pm 10\%$ of the setting value.

- The measurement value is within $\pm 10\%$ of the setting value.
Proceed to *Step 5*).
- The measurement value is not within $\pm 10\%$ of the setting value.

Change the parameter for “Motor Pulse” and press **Enter**.


Repeat from *Step 3*) to *Step 4*) until the value is within the criteria.

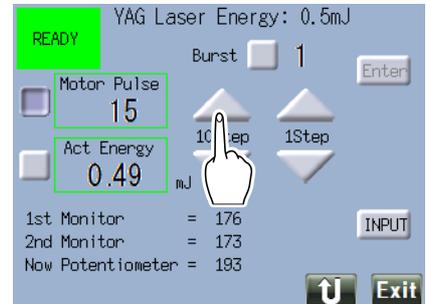

5) Press “Act Energy”.



6) Enter the measurement value after the calibration.

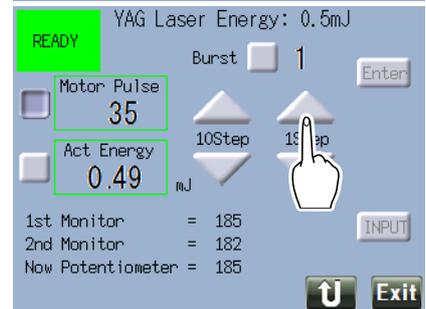

7) Press [Input] to save the parameters for “Motor Pulse” and “Act Energy”.



Knowledge

- [Input] can be pressed after the laser is emitted.



8) Change “Burst” to 2.

Repeat the treatment beam emission and change the parameter for “Motor Pulse” until the measurement value is within $\pm 10\%$ of the setting value. Enter the measurement value after the calibration (the sum of the 2 pulses) to the “Act Energy”.


When the measurement value is within the criteria, press [Input] to save the parameter.



9) Change "Burst" to 3.

Repeat the treatment beam emission and change the parameter for "Motor Pulse" until the measurement value is within $\pm 10\%$ of the setting value. Enter the measurement value after the calibration (the sum of the 3 pulses) to the "Act Energy".

When the measurement value is within the criteria, press [Input] to save the parameter.

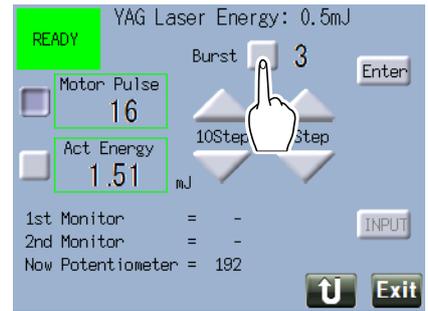


10) Press .

7 Change the parameters for 1.0 mJ to MAX energy values.

Change the parameters in the same manner as with 0.5 mJ energy value, starting from the smallest energy value.

Concerning "MAX", perform calibration for the following combinations.

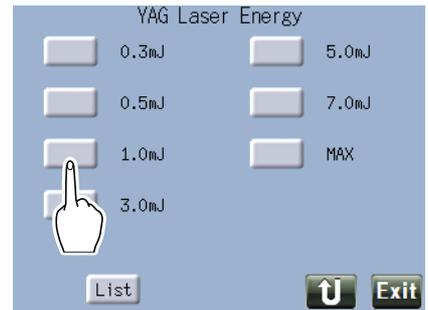

Mode	ENERGY	BURST
YAG mode	8.0 mJ	3
	9.0 mJ	2
	10.0 mJ	1

8 Record the measurement values.

 "5.5 Laser Beam Output Calibration Record" (page 98)

9 Press .



Knowledge

- In SLT mode, calibrate the treatment beam for 0.3 mJ, 0.5 mJ, 1.0 mJ, 2.0 mJ, and 3.0 mJ energy values.
In the same manner as YAG mode, change the parameters in the order starting from the smallest energy value.
In SLT mode, the Burst button is not displayed.

4.8 Measuring and Calibrating Aiming Beam Power Output

- If the power output of the aiming beam is out of tolerance, calibration is required.
A power meter is necessary for power output measurement and calibration.


CAUTION

- Only service personnel trained by NIDEK or NIDEK distributor are allowed to calibrate the laser power output.
- Follow the procedure at least once a year.

4.8.1 Measuring aiming beam power output


1 Mount a power meter.

- 1) Secure the detector of the power meter to the chinrest.
Align the height of the center of the detector's receiving surface with the eye level marker.
- 2) Connect the detector cable to the power meter.

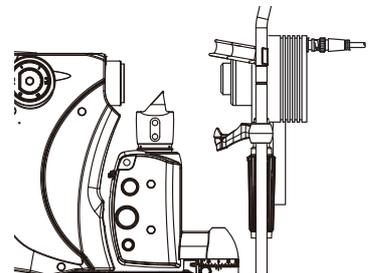


2 Activate the device.

 "3.2.1 Activating the device" (page 41)

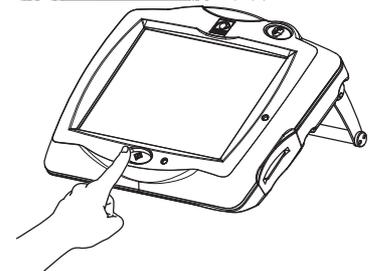


3 Align the aiming beam with the receiving surface of the detector.

- 1) Turn off the illumination light by turning the illumination control.
- 2) Project the aiming beam on the receiving surface of the detector.

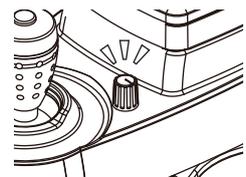

- 3) Manipulate the joystick to align the aiming beam with a position where the aiming beam spot is within the receiving surface.
It is not necessary to adjust the focus.

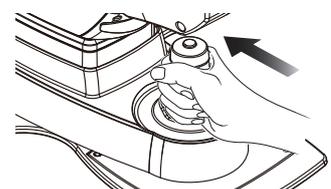

4 Measure the aiming beam for the combinations in the table.

Mode	Setting	Tolerance
YAG mode	1	$0.5 \pm 0.3 \mu\text{W}$
	15	$25 \pm 5 \mu\text{W}$
SLT mode	1	0.04 mW or less
	15	$0.3 \pm 0.1 \text{ mW}$

5 Record the measurement values.

 "5.5 Laser Beam Output Calibration Record" (page 98)

6 When the measurement values are out of the tolerance, the power output needs to be calibrated.


 "4.8.2 Calibrating aiming beam power output" (page 86)

4.8.2 Calibrating aiming beam power output

1 Activate CALIBRATION mode.

↳ “3.4 Changing Device Settings” (page 50)

Only administrators qualified by NIDEK are informed of the password.

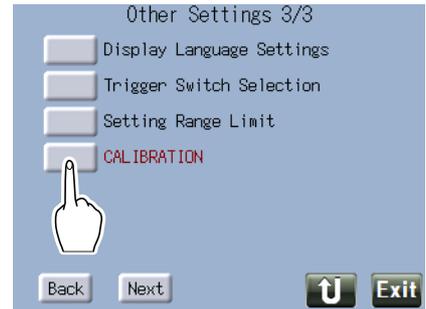

2 Press “Aiming Power”.



3 Change the parameters.

1) Press a button to change the corresponding parameter.

Minimum power output: AIMING OFFSET

Maximum power output: AIMING GAIN



Knowledge

- The tolerance of the measurement values compared to the setting values is described below.
- YAG mode
Maximum: $25 \pm 5 \mu\text{W}$
Minimum: $0.5 \pm 0.3 \mu\text{W}$
- SLT mode
Maximum: $0.3 \pm 0.1 \text{ mW}$
Minimum: 0.04 mW or less



2) Press **Enter** .

3) Press **Up** .

4 Record the measurement values.

↳ “5.5 Laser Beam Output Calibration Record” (page 98)

5 Press **Exit** .



SPECIFICATIONS AND TECHNICAL INFORMATION

5.1 Specifications

Device specifications		
YAG treatment beam	• Wavelength	1,064 nm
	• Energy output	BURST 1: 0.3 to 10.0 mJ ±20% BURST 2: 0.3 to 9.0 mJ ±20% BURST 3: 0.3 to 8.0 mJ ±20%
	• Energy output stability	The energy output meets the specification 90% or greater (out of continuous 20 shots).
	• Number of emitted pulses	Single mode: one shot / trigger Burst mode: two or three shots / trigger
	• Spot size	8 µm
	• Pulse width	3 ns ±20%
	• Pulse interval	60 to 200 µs
	• Continuous irradiation / repetition time	Single mode: 3 Hz Burst mode: 1.5 Hz
	• Cone angle	16° ±2°
	• Focus shift	±500 µm
	• Nominal Ocular Hazard Distance (NOHD)	Single mode: 2.65 m Burst mode: 1.71 m
YAG aiming beam	• Wavelength	635 nm
	• Power output	OFF, 0.5 to 25 µW Max. 25 ±5 µW, Min. 0.5 ±0.3 µW
Protection of operator's eye (YAG)	• Optical density (O.D.)	5 or more

Device specifications		
SLT treatment beam (YC-200 S plus)	• Wavelength	532 nm
	• Energy output	0.3 to 3.0 mJ \pm 20%
	• Energy output stability	The energy output meets the specification 90% or greater (out of continuous 120 shots).
	• Spot size	400 μ m
	• Pulse width	3 ns \pm 20%
	• Continuous irradiation / repetition time	Single mode: 3 Hz
	• Cone angle	5.5° \pm 1°
	• Nominal Ocular Hazard Distance (NOHD)	13.7 m
SLT aiming beam (YC-200 S plus)	• Wavelength	635 nm
	• Power output	OFF, 0.04 to 0.3 mW Max. 0.3 \pm 0.1 mW, Min. 0.04 mW or less
	• Cone angle	7.5° or less
	• Spot size	400 μ m
Protection of operator's eye (SLT) (YC-200 S plus)	• Optical density (O.D.)	5 or more
Microscope optical system	• Optical axis	Galilean type
	• Distance between right and left optical axes	22 mm
	• Objective lens focal length	f = 125 mm
	• Total magnification	5x, 8x, 12.5x, 20x, 32x
	• Real field of view	\varnothing 40.7 mm, 25.7 mm, 16.1 mm, 10.1 mm, 6.4 mm
Binocular microscope	• Optical system	Convergence angle 6°
	• Distance between right and left optical axes	22 mm
	• Minimum interpupillary adjustment range	50 to 78 mm
	• Eye relief	20 mm or more
	• Eyepiece magnification	12.5x
	• Lens aperture	23 mm
	• Diopter adjustment range	\pm 8 D

Device specifications		
Slit illumination	• Illumination light source	White LED
	• Slit length	14 mm, 9 mm, 5 mm, 3 mm, 1 mm, 0.2 mm
	• Slit width	0 to 14 mm
	• Slit rotation	360°
	• Illuminance	OFF, minimum illuminance to maximum illuminance (continuously variable) Max. 208,000 lx to 160,000 lx * Illuminance at 0° when the illumination tower equipped with a tilting function is used
	• Filter	Blue, red-free, 15% ND
Illumination tower (SLT) Illumination tower equipped with a tilting function (YC-200 S plus)	• Illumination incidence angle	0 to 20° (downward)
Illumination tower (YAG) Illumination tower with the base fixed	• Illumination incidence angle	18° (downward)
Head rest	• Eye level height	395 mm
	• Distance between supports	220 mm
	• Chinrest up and down movement	80 mm (manual)
Fixation lamp	• Light source	LED (green)
Main body horizontal movement	• Forward and backward	80 mm
	• Right and left	100 mm
	• Fine movement distance	±5 mm
Main body up and down movement	• Up and down	30 mm (motorized)
Power supply	• Voltage	AC 100 to 240 V ±10%
	• Frequency	50/60 Hz
	• Power consumption	A maximum of 100 VA
Dimensions and mass	• Dimensions	346 mm (W) × 422 mm (D) × 577 mm (H)
	• Mass	17.0 kg (YC-200) 18.0 kg (YC-200 S plus)

Device specifications		
Environmental conditions (during use)	• Temperature	15 to 30°C (59 to 86°F)
	• Humidity	30 to 90% (non-condensing)
	• Atmospheric pressure	800 to 1,060 hPa
	• Installation location	Indoors On a level, stable location in a darkened room or laser treatment room with laser protective enclosure
	• Other	No harmful dust or smoke
Environmental conditions (during transport [packed condition] and storage)	• Temperature	-10 to 55°C (14 to 131°F)
	• Humidity	10 to 95% (non-condensing)
	• Atmospheric pressure	500 to 1,060 hPa
Expected service life	• 7 years (defined by manufacturer)* Proper maintenance is necessary.	
Equipment to be used in combination	• NIDEK ophthalmic laser photocoagulator Green Laser Photocoagulator GYC-500	
Classifications	<ul style="list-style-type: none"> • Laser classification (IEC 60825-1): Class 3B Laser product that are normally hazardous when intrabeam ocular exposure occurs within the NOHD including accidental short time exposure. It is necessary to avoid exposure to laser. Viewing diffuse reflections is normally safe. • Protection against electrical shock: Class I ME equipment • Protection against electrical shock (applied parts): Type B applied part • Protection against harmful ingress of water or particulate matter: IPX0 (main body), IPX8 (optional foot switch) • Method(s) of sterilization: ME equipment that does not contain any part that needs sterilization. • Suitability for use in an oxygen rich environment: ME equipment that is not intended for use in an oxygen rich environment • Mode of operation: Continuous operation • Conformity to the standard for slit-lamp microscope: In compliance with ISO 10939:2007 	

Device specifications		
Accessories		
Standard accessories	Main body, head rest, control box, connector box, remote connector plug, keycard, power cord, chinrest paper, head belt, grips, focusing rod, attachment plates, arm rest, cap, dust cover, DANGER label, operator's manual	
Optional accessories	Foot switch, optical table, illumination tower (tilting / with the base fixed / split mirror illumination), stand for CB, safety goggles	
Optional optical table specifications		
Operation	• Elevation range	653 to 851 mm (Table top height)
	• Elevation time	12 s (60 Hz)
Power supply	• Voltage	For 100 V regions: AC 120 V \pm 10% For 200 V regions: AC 230 V \pm 10%
	• Frequency	50/60 Hz
	• Power consumption	1000VA
Dimensions and mass	• Dimensions	742 mm (W) \times 463 mm (D) \times 653 to 851 mm (H)
	• Mass	25 kg
Environmental conditions	• The same as those of the main body specifications	

- If you need information to ensure cybersecurity, contact Nidek or your authorized distributor.


5.2 Light Hazard

- Provision of information on the avoidance of light hazard from the optical device is required in ISO 15004-2:2007 "Ophthalmic instruments - Fundamental requirements and test methods -".




CAUTION

- The light emitted from this device is potentially hazardous. The longer the duration of exposure is, the greater the risk of ocular damage becomes.
Exposure to light from this device when operated at maximum intensity will exceed the safety guideline at the following duration (except for the treatment beam).
During treatment of posterior capsule (YC-200): 303 seconds
During treatment of anterior chamber angle (YC-200 S plus): 347 seconds

- Graphs of relative spectral power output

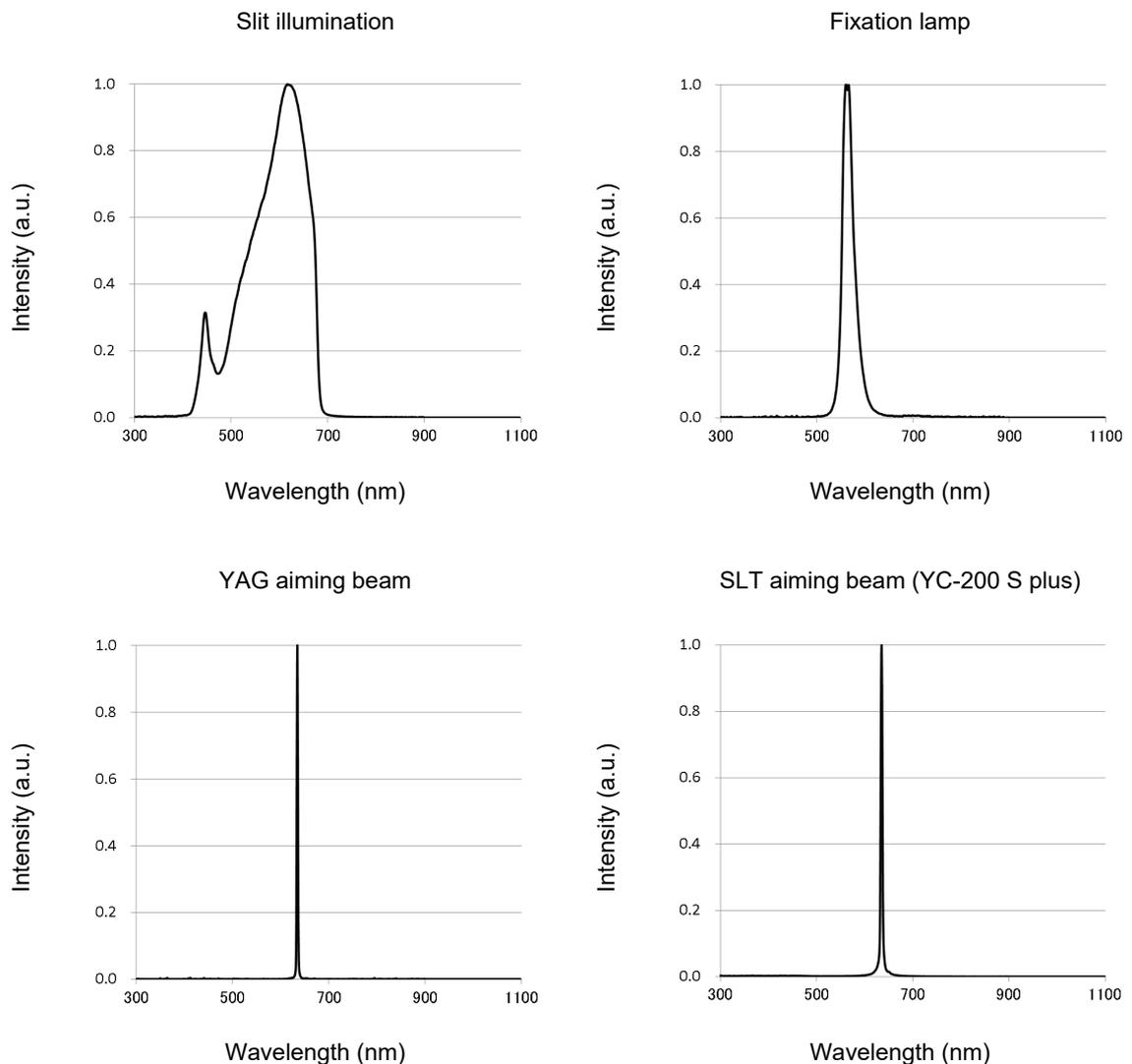

● Maximum intensity and ratio to maximum intensity

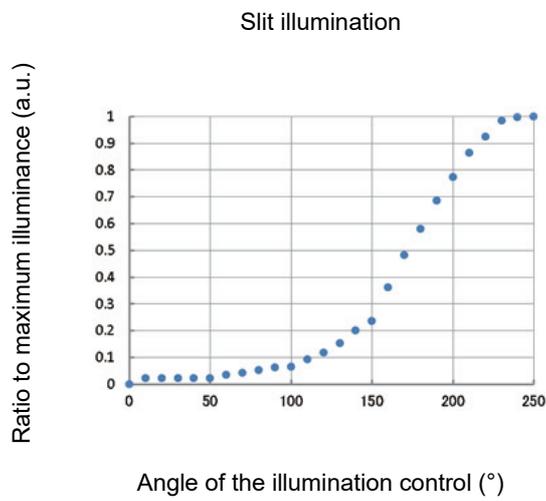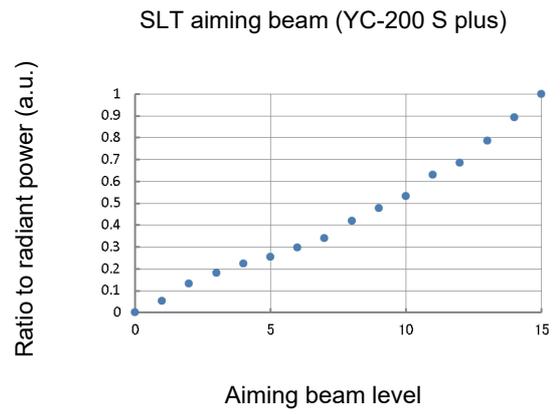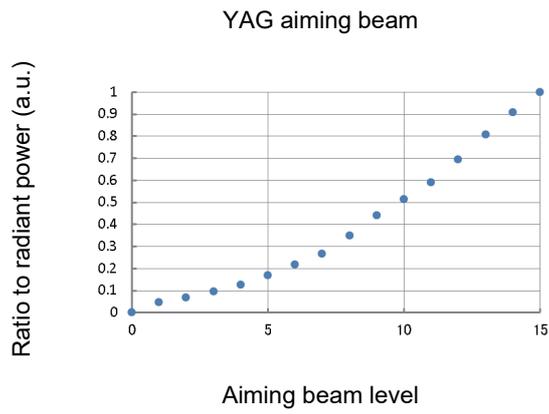

5.3 EMC (Electromagnetic Compatibility)

- The device is suitable for use in hospitals except for near active HF surgical equipment and RF shielded rooms with an ME system for magnetic resonance imaging, where the intensity of electromagnetic disturbances is high, electrophysiology laboratories, or areas where short-wave therapy equipment is used.

WARNING

- Do not use the device near, on, or under other electronic equipment or electromagnetic disturbance sources. Otherwise, it could result in improper operation. If such use is necessary, the device and the other equipment should be observed to verify that they are operating normally.
- Use of accessories, cables other than those specified or provided by the manufacturer of this equipment could result in increased electromagnetic emissions or decreased electromagnetic immunity of this equipment and result improper operation.
- Portable RF communications equipment (including peripherals such as antenna cables and external antennas) or electromagnetic disturbance sources as shown below should be used no closer than 30 cm (12 inches) to any part of the device, including the specified or provided cables. Otherwise, degradation of the performance of this equipment could result.

The following are examples of electromagnetic disturbance sources:

- Induction cooking appliance and ovens
- RFID readers
- Electronic article surveillance (EAS) systems
- Sponge detection systems
- Equipment used for position detection (e.g. in catheter labs)
- Wireless power transfer charging systems for electrical vehicles

◆ Specified cable

Part name	Cable shielded	Ferrite core	Length (m)
Communications cable for CB	Yes	No	1.1
Main body connection cable	Yes	No	0.55
Fixation lamp connection cable	Yes	No	0.7
YC connection cable	Yes	No	1.3
GYC connection cable	Yes	Yes	2.9
Fiber optic cable (COMBO)	No	No	3.0
Foot switch cable	Yes	No	2.9
Power cord	No	No	3.0

◆ Essential performance

- Laser irradiation function

◆ Compliance for Emission Standard

Phenomenon	Product family standard	Compliance
Conducted and radiated RF emissions	CISPR 11	Group 1 Class B
Harmonic distortion	IEC 61000-3-2	*1
Voltage fluctuations and flicker	IEC 61000-3-3	*2

* 1 For the regions where the rated voltage is 220 V to 240 V, this device complies with this standard.

* 2 For the regions where the rated voltage (line to neutral) is 220 V to 250 V, this device complies with this standard.

◆ Test specifications for enclosure port immunity to RF wireless communications equipment

Test frequency (MHz)	Band (MHz)	Service	Modulation	Immunity test level (V/m)
385	380 to 390	TETRA 400	Pulse modulation 18 Hz	27
450	430 to 470	GMRS 460, FRS 460	FM ±5 kHz deviation 1 kHz sine	28
710	704 to 787	LTE Band 13, 17	Pulse modulation 217 Hz	9
745				
780				
810	800 to 960	GSM 800/900, TETRA 800, iDEN 820, CDMA 850, LTE Band 5	Pulse modulation 18 Hz	28
870				
930				
1720	1700 to 1990	GSM 1800; CDMA 1900; GSM 1900; DECT; LTE Band 1, 3, 4, 25; UMTS	Pulse modulation 217 Hz	28
1845				
1970				
2450	2400 to 2570	Bluetooth WLAN 802.11 b/g/n RFID 2450 LTE Band 7	Pulse modulation 217 Hz	28
5240	5100 to 5800	WLAN 802.11 a/n	Pulse modulation 217 Hz	9
5500				
5785				

◆ Compliance for Immunity Standard

Phenomenon	Basic EMC standard	Immunity test levels
Electrostatic discharge	IEC 61000-4-2	±8 kV contact ±2 kV, ±4 kV, ±8 kV, ±15 kV air
Radiated RF electromagnetic field	IEC 61000-4-3	3 V/m 80 MHz - 2.7 GHz 80% AM at 1 kHz
Proximity fields from RF wireless communications equipment	IEC 61000-4-3	See "Test specifications for enclosure port immunity to RF wireless communications equipment".
Electrical fast transients / bursts	IEC 61000-4-4	Input power port ±2 kV 100 kHz repetition frequency
		Signal input/output parts port ±1 kV 100 kHz repetition frequency
Surges Line-to-line	IEC 61000-4-5	Input power port ±0.5 kV, ±1 kV
Surges Line-to-ground		Input power port ±0.5 kV, ±1 kV, ±2 kV Signal input/output parts port ±2 kV
Conducted disturbances induced by RF fields	IEC 61000-4-6	3 V 0.15 MHz – 80 MHz 6 V in ISM bands between 0.15 MHz and 80 MHz 80% AM at 1 kHz
Rated power frequency magnetic fields	IEC 61000-4-8	30 A/m 50 Hz or 60 Hz
Voltage dips	IEC 61000-4-11	0% U _T ; 0.5 cycle At 0°, 45°, 90°, 135°, 180°, 225°, 270°, and 315°
		0% U _T ; 1 cycle and 70% U _T ; 25/30 cycles Single phase: at 0°
Voltage interruptions	IEC 61000-4-11	0% U _T ; 250/300 cycles

5.5 Laser Beam Output Calibration Record

■ Print it if necessary.

● YAG treatment beam

Setting	Measurement value (mJ)		
	Single (Burst 1)	Burst 2	Burst 3
0.3 mJ			
0.5 mJ			
1.0 mJ			
3.0 mJ			
5.0 mJ			
7.0 mJ			
MAX			

● SLT treatment beam

Setting	Measurement value (mJ)
0.3 mJ	
0.5 mJ	
1.0 mJ	
2.0 mJ	
3.0 mJ	

● YAG aiming beam

Setting	Tolerance	Measurement value (μW)
1	$0.5 \pm 0.3 \mu\text{W}$	
15	$25 \pm 5 \mu\text{W}$	

● SLT aiming beam

Setting	Tolerance	Measurement value (mW)
1	0.04 mW or less	
15	0.3 ±0.1 mW	

5.6 Pre-use Checklist

- Print it if necessary.

Item	Date and person responsible for check
The lenses and mirrors are clean.	
The forehead rest, chinrest, grips, and head belt are clean.	
The power cord is connected to a grounded power outlet that meets the power requirements.	
The power indicator on the control box illuminates.	
Inserting the keycard and pressing the start button turns on the screen of the control box.	
YAG mode	
The AIMING value increases and decreases.	
AIMING rotates.	
The ENERGY value increases and decreases. After the change, the test fire is performed.	
The FOCUS SHIFT value increases and decreases.	
The BURST value changes. After the change, the test fire is performed.	
After AIMING lights up, "READY" and "STANDBY" can be switched.	
After the device enters READY mode, pressing the hand switch emits the treatment beam.	
The aiming beam and treatment beam are coaxial.  (page 73)	
Pressing COUNT clears the value.	
The device is turned off when the emergency laser stop button is pressed.	
SLT mode (YC-200 S plus)	
Pressing  toggles YAG mode and SLT mode.	
The AIMING value increases and decreases.	
The ENERGY value increases and decreases. After the change, the test fire is performed.	
SLT-NAVI can be set.	
After AIMING lights up, "READY" and "STANDBY" can be switched.	
After the device enters READY mode, pressing the hand switch emits the treatment beam.	
The aiming beam and treatment beam are coaxial.  (page 75)	
Pressing COUNT clears the value.	
The device is turned off when the emergency laser stop button is pressed.	

5.7 Glossary

Term	Details
Device	The main laser device of YAG laser system Two types are available: YC-200 and YC-200 S plus
YC-200	The device equipped with YAG mode function only.
YC-200 S plus	The device equipped with YAG and SLT mode functions.
CB	Control box
Connector box	A unit that contains the power supply and connectors of the device.
YAG mode	An operation mode in which treatment using the YAG treatment beam (wavelength: 1,064 nm) is available. It is used mainly for posterior capsulotomy and iridotomy.
SLT mode	An operation mode in which treatment using the SLT treatment beam (wavelength: 532 nm) is available. It is used for selective laser trabeculoplasty.
Laser beam	Treatment beam and aiming beam
YAG laser YAG treatment beam	Laser beam for treatment of the affected tissue in YAG mode. A 1,064 nm pulsed laser beam is used.
SLT laser SLT treatment beam	Laser beam for treatment of the affected tissue in SLT mode. A 532 nm pulsed laser beam (green) is used.
Laser energy output	Laser energy emitted from the laser aperture of the device. (Unit: mJ)
YAG aiming beam	Laser beam that indicates the position to which the treatment beam is to be emitted in YAG mode. The two-beam system (separating the YAG aiming beam into two beams) is used. The focus position is determined according to the alignment of the beams.
SLT aiming beam	Laser beam that indicates the position to which the treatment beam is to be emitted in SLT mode. The parfocal optical system is used. The focus position is determined according to the projection status of the beams.
Test fire	Checking that there are no abnormalities with the laser by emitting the laser within the device before the laser is actually emitted toward the patient's eye
Spot size	Diameter of the laser beam spot
Single mode	One shot of treatment beam is emitted each time the trigger switch is pressed in this mode.
Burst mode (YAG mode)	Two or three shots of treatment beam are emitted each time the trigger switch is pressed in this mode.
Focus shift (YAG mode)	To shift the focal points (focus) of the YAG treatment beam and YAG aiming beam. It is mainly used for the treatment of secondary cataract in which the focal point of the YAG treatment beam is shifted to the posterior chamber (POST.) side from the aiming beam to prevent inclusion of pit in the intraocular lens.
Aiming beam rotation (YAG mode)	To rotate the spot order of the YAG aiming beam in order to set the spots to the desired positions. The function of continuous rotation for the specified durations is also available.
SLT-NAVI (SLT mode)	The function that assists the operator in surgery by specifying the laser emission positions and orders before the treatment Enter the emission area, emission start position, and emission direction.

Term	Details
READY	The treatment beam can be emitted. Pressing the status button of the control box toggles between "STANDBY" and "READY".
STANDBY	The treatment beam cannot be emitted. Pressing the status button of the control box toggles between "STANDBY" and "READY".
Protective filter	The optical filter to protect the operator's eye from the reflected light of the treatment beam Inserted and fixed within the device. YC-200: Only YAG treatment beam (1,064 nm) is blocked. YC-200 S plus: YAG treatment beam (1,064 nm) and SLT treatment beam (532 nm) are blocked.
Grips	Patients hold them to keep themselves in a stable posture during surgery.
Chinrest elevation control	Used to adjust the height of the patient's chin.
Chinrest	Used to rest the patient's chin.
Fixation lamp	Used to steady the patient's visual axis by having the patient focus on the lamp. The position of the fixation lamp can be adjusted by the flexible arm.
Forehead rest	Used to rest the patient's forehead to keep them in a stable posture.
Eye level marker	Markers to adjust the height of the patient's eyes.
Base unit fastening knob	Used to fasten the device so that the device may not move horizontally.
Illumination unit arm fastening knob	Used to secure the illumination unit arm.
Microscope arm fastening knob	Used to secure the microscope arm.
Eyepieces	Lenses attached to the part closest to the operator's eyes.
Magnification changer	Used to select the total magnification of the microscope.
Slit rotation control	Used to rotate the illumination light.
Slit width control	Used to continuously adjust the width of the illumination light.
Slit length control	Used to select the aperture diameter (slit length) of the illumination light.
Filter changer	A disk to select a filter to be inserted into the observation optical path
Joystick	Used to adjust the position of the observation unit. By manipulating the joystick right, left, forward, and backward, the unit moves right, left, forward, and backward. By tilting the joystick, the unit slightly moves to the direction of tilt.
Hand switch	Used to emit the treatment beam. Select either the hand switch or the optional foot switch as a trigger switch.
S-Switch	An auxiliary switch built-in the joystick. Three functions can be assigned as desired.
Illumination control	Used to adjust the intensity of the illumination light.
Power indicator	Lights up when the master switch is turned on while power is supplied.
Emergency laser stop button	An emergency button to stop the power and operation of the device After the stop, pressing the start button restarts the device.
Start button	A button to activate and stop the device
Keycard	A memory card equipped with a function as a key necessary to activate the device The configuration information of the device and error log information are also recorded.

Term	Details
Access indicator	An indicator that indicates whether the keycard is inserted and whether data is read or written in the keycard. It lights up when the keycard is inserted and blinks when data is read or written in the keycard.
Touch screen	A screen to specify and display conditions
Limit mode (YAG mode)	A function to limit the setting ranges of Energy, Focus Shift, and Burst to use the device more safely
REMOTE connector	A connector designed to turn off the device when an external switch connected by a signal line is turned off. Connect this connector to an external switch such as the door switch of the operating room so that the device is stopped when an outsider enters the operating room. In addition, this function allows a surgeon other than the one conducting the surgery to stop the device by turning off the external switch connected to the connector based on their judgment in case of a danger during surgery.
Nominal Ocular Hazard Distance (NOHD)	The Nominal Ocular Hazard Distance (NOHD) is the distance along the axis of the unobstructed beam from the laser aperture where the irradiance falls below the applicable exposure limit.
Parfocal optical system (YC-200 S plus)	An optical system in which the image of an object surface is formed on the target surface. The SLT aiming beam is emitted from the fiber tip (the object surface) in the parfocal optical system so that it appears as a sharply-edged spot on the target surface.
Expected service life	A period of time beyond which the reliability and safety of the device cannot be guaranteed even under normal use and regular maintenance that involves repeated replacement of maintenance parts and consumable parts, repair, and overhaul.



INDEX

A

Access indicator	28
AIMING	.31, 34

B

Base unit fastening knob	25
BURST	31
Burst mode	101

C

CB	101
Chinrest	23
Chinrest elevation control	23
Cleaning	71
Connector box	26
COUNT	.30, 33

E

Emergency laser stop button	27
ENERGY	.32, 34
Eye level marker	23
Eyepieces	102

F

Filter changer	25
Fixation lamp	23
FOCUS SHIFT	32
Focus shift	101
Forehead rest	23

G

Grips	23
-------	----

H

Hand switch	24
-------------	----

I

Illumination control	24
Illumination unit arm fastening knob	25
Information display field	.32, 34

J

Joystick	23
----------	----

K

Keycard	21
Keycard slot	28

L

Laser beam	101
Laser energy output	101
Limit mode	55

M

Magnification changer	24
Master switch	26
MENU button	.32, 35
Microscope	23
Microscope arm fastening knob	25
Mode change button	.32, 35

P

Power indicator	27
Power inlet	26
Prism head	24
Protective filter	102

R

READY	102
REMOTE connector	26

S

Single mode	101
Slit length control	25
Slit rotation control	25
Slit width control	25
SLT aiming beam	101
SLT laser	101
SLT mode	17, 33, 101
SLT treatment beam	101
SLT-NAVI	.34, 48
Spot size	.87, 88
S-Switch	25
STANDBY	102
Start button	27
Status button	.30, 33
Status indicator	27

T

Test fire45, 101
Σ ENERGY31, 34
Touch screen	27

Y

YAG aiming beam	101
YAG laser	101
YAG mode	17, 30, 36, 101
YAG treatment beam	101
YC-200	17
YC-200 S plus	17